

VAMAS Inter-laboratory study on intensity calibration for XPS instruments using low-density polyethylene

Running title: VAMAS study on intensity calibration for XPS instruments using LDPE

Running Authors: Reed et al.

Benjamin P. Reed^{*1}, David J. H. Cant¹, Steve J. Spencer¹, Abraham Jorge Carmona-Carmona², Adam Bushell³, Alberto Herrera-Gómez², Akira Kurokawa⁴, Andreas Thissen⁵, Andrew G. Thomas⁶, Andrew J. Britton⁷, Andrzej Bernasik⁸, Anne Fuchs⁹, Arthur P. Baddorf¹⁰, Bernd Bock¹¹, Bill Theilacker¹², Bin Cheng¹³, David G. Castner¹⁴, David J. Morgan¹⁵, David Valley¹⁶, Elizabeth A. Willneff¹⁷, Emily F. Smith¹⁸, Emmanuel Nolot¹⁹, Fangyan Xie²⁰, Gilad Zorn²¹, Graham C. Smith²², Hideyuki Yasufuku²³, Jeffery Fenton²⁴, Jian Chen²⁰, Jonathan D. P. Counsell²⁵, Jörg Radnik²⁶, Karen J. Gaskell²⁷, Kateryna Artyushkova¹⁶, Li Yang²⁸, Lulu Zhang⁴, Makiho Eguchi²⁹, Marc Walker³⁰, Mariusz Hajdyla⁸, Mateusz M. Marzec⁸, Matthew R. Linford³¹, Naoyoshi Kubota²⁹, Orlando Cortazar-Martínez², Paul Dietrich⁵, Riki Satoh²⁹, Sven L. M. Schroeder⁷, Tahereh G. Avval³¹, Takaharu Nagatomi³², Vincent Fernandez³³, Wayne Lake³⁴, Yasushi Azuma⁴, Yusuke Yoshikawa^{35,36}, and Alexander G. Shard¹

*Corresponding author: Tel: +44 (0)20 8943 6952; E-mail: benjamin.reed@npl.co.uk

¹National Physical Laboratory, Hampton Road, Teddington, TW11 0LW, UK.

²CINVESTAV-Unidad Queretaro, Queretaro 76230, Mexico.

³Thermo Fisher Scientific (Surface Analysis), East Grinstead, RH19 1XZ, UK.

⁴National Metrology Institute of Japan (NMIJ), National Institute of Advanced Industrial Science and Technology (AIST), 1-1-1 Higashi, Tsukuba, Ibaraki 305-8565, Japan.

⁵SPECS Surface Nano Analysis GmbH, Voltastraße 5, 13355 Berlin, Germany.

⁶Department of Materials, Photon Science Institute and Sir Henry Royce Institute, Alan Turing Building, University of Manchester, Oxford Road, Manchester, M13 9PL, UK.

⁷Versatile X-ray Spectroscopy Facility, School of Chemical and Process Engineering, University of Leeds, Leeds, LS2 9JT, UK.

⁸Academic Centre for Materials and Nanotechnology, AGH University of Science and Technology, 30-059 Kraków, Poland.

⁹Robert Bosch GmbH, Robert-Bosch-Campus, 71272 Renningen, Germany.

¹⁰Center for Nanophase Materials Sciences, Oak Ridge National Laboratory, 1 Bethel Valley Rd, Oak Ridge, TN 37830, USA.

¹¹Tascon GmbH, Mendelstr. 17, D-48149 Münster, Germany.

¹²Medtronic, 710 Medtronic Parkway, LT240, Fridley, MN 55432, USA.

- ¹³Analysis and Testing Center, Beijing University of Chemical Technology, Beijing 100029, P. R. China.
- ¹⁴National ESCA and Surface Analysis Center for Biomedical Problems, Department of Bioengineering and Chemical Engineering, University of Washington, Seattle, Washington 98195, USA.
- ¹⁵Cardiff Catalysis Institute, School of Chemistry, Cardiff University, Main Building, Cardiff CF10 3AT, UK.
- ¹⁶Physical Electronics Inc., East Chanhassen, Minnesota 55317, USA.
- ¹⁷Versatile X-ray Spectroscopy Facility, School of Design, University of Leeds, Leeds, LS2 9JT, UK.
- ¹⁸Nanoscale and Microscale Research Centre, University of Nottingham, Nottingham, NG7 2QB, UK.
- ¹⁹CEA-LETI, 17 rue des Martyrs, 38054 Grenoble, France.
- ²⁰Instrumental Analysis & Research Center, Sun Yat-sen University, Guangzhou 510275, P. R. China.
- ²¹GE Research, 1 Research Circle, K1 1D7A, Niskayuna, NY 12309, USA.
- ²²Faculty of Science and Engineering, University of Chester, Thornton Science Park, Chester, CH2 4NU, UK.
- ²³Materials Analysis Station, National Institute for Materials Science (NIMS), 1-2-1 Sengen, Tsukuba, Ibaraki 305-0044, Japan.
- ²⁴Medtronic, 6700 Shingle Creek Parkway, Brooklyn Center, MN 55430, USA.
- ²⁵Kratos Analytical Ltd., Wharfside, Trafford Wharf Road, Manchester M17 1GP, UK.
- ²⁶Bundesanstalt für Materialforschung und -prüfung (BAM), Unter den Eichen 44-46, 12203 Berlin, Germany.
- ²⁷Department of Chemistry and Biochemistry, University of Maryland, College Park, MD 20742, USA.
- ²⁸Department of Chemistry, Xi'an Jiaotong-Liverpool University, Suzhou Dushu Lake Higher Education Town, P. R. China.
- ²⁹Surface Analysis Dep. Materials Characterization Div., Futtsu Unit. NIPPON STEEL TECHNOLOGY Co.,Ltd., Japan.
- ³⁰Department of Physics, University of Warwick, Coventry, West Midlands, CV4 7AL, UK.
- ³¹Department of Chemistry and Biochemistry, Brigham Young University, Provo, UT 84604, USA.
- ³²Platform Laboratory for Science and Technology, Asahi Kasei Corporation, 2-1 Samejima, Fuji, Shizuoka 416-8501, Japan.
- ³³Institut des matériaux Jean Rouxel, 44300 Nantes, France.
- ³⁴Atomic Weapons Establishment (AWE), Aldermaston, Reading, Berkshire, RG7 4PR, UK.
- ³⁵Material Analysis Department, Yazaki Research and Technology Center, YAZAKI Corporation, Susono City, Japan.
- ³⁶Technical & Project Support Dept., Nidec Sankyo CMI Corporation, 46-1, Senpuku, Susono-city, Shizuoka 460-1116, Japan (current address).

We report the results of a VAMAS (Versailles Project on Advanced Materials and Standards) inter-laboratory study on the intensity scale calibration of X-ray photoelectron spectrometers using low-density polyethylene (LDPE) as an alternative material to gold, silver, and copper. An improved set of LDPE reference spectra, corrected for different instrument geometries using a quartz-monochromated Al K α X-ray source, was developed using data provided by participants in this study. Using these new reference spectra, a transmission function was calculated for each dataset that participants provided. When compared to a similar calibration procedure using the NPL reference spectra for gold, the LDPE intensity calibration method achieves an absolute offset of $\sim 3.0\%$ and a systematic deviation of $\pm 6.5\%$ on average across all participants. For spectra recorded at high pass energies (≥ 90 eV), values of absolute offset and systematic deviation are $\sim 5.8\%$ and $\pm 5.7\%$ respectively, whereas for spectra collected at lower pass energies (< 90 eV), values of absolute offset and systematic deviation are $\sim 4.9\%$ and $\pm 8.8\%$ respectively; low pass energy spectra perform worse than the global average, in terms of systematic deviations, due to diminished count rates and signal-to-noise. Differences in absolute offset are attributed to the surface roughness of the LDPE induced by sample preparation. We further assess the usability of LDPE as a secondary reference material and comment on its performance in the presence of issues such as: variable dark noise, X-ray warm up times, inaccuracy at low count rates, and underlying spectrometer problems. In response to participant feedback and the results of the study, we provide an updated LDPE intensity calibration protocol to address the issues highlighted in the inter-laboratory study. We also comment on the lack of implementation of a consistent and traceable intensity calibration method across the community of XPS users, and therefore propose a route to achieving this with the assistance of instrument manufacturers, metrology laboratories, and experts leading to an international standard for XPS intensity scale calibration.

I. INTRODUCTION

X-ray photoelectron spectroscopy (XPS) is a highly surface sensitive technique that can determine the elemental composition and chemical environment of the top few nanometres (~ 10 nm) of a material's surface. XPS has become synonymous with surface analysis, with many industries and academic

institutes regularly using quantitative XPS in their research, especially as it becomes more accessible through commercial XPS instruments.¹ Quantitative XPS analysis relies on the extraction of accurate peak intensities (i.e. their area) from measured spectra in order to ascertain the concentration of specific elements and chemical species in the sample. The intensity, I_n , of a photoelectron peak for a given element, n , is affected by both physical factors and experimental factors such that, for a homogeneous material,

$$I_n \sim N_n J \lambda_E F \frac{\sigma}{4\pi} T \quad (1)$$

where N_n is the atomic concentration of element n , J is the X-ray flux incident on the sample, λ_E is the energy-dependent inelastic mean free path,² F is a term related to the photoelectron angular emission distribution and analyser geometry,³ σ is the core level photoionization cross-section,^{4,5} and T is the transmission efficiency of photoelectrons through the analyser as a function of kinetic energy, commonly called the ‘transmission function’. Except for N_n (assuming it is unknown) and T , all these factors can be estimated through theoretical calculations, found in look-up tables, or factored out when comparing relative intensities acquired using the same experimental parameters. As with all experimental apparatus, XPS instruments require regular calibration. The intensity scale of an XPS instrument is calibrated by determining the T of the analyser, which is affected by the detector efficiency and factors such as the overlap between the X-ray spot and the collection area of the analyser for a given sample positioning and orientation (usually, the position of the sample is chosen to maximize counts; the calibration is usually done with the sample facing the analyser). T is also affected by analyser settings such as the pass energy, defined sample analysis area, slit/iris settings, voltage across the micro-channel plate, and the operating mode of the electromagnetic lens stack, to name a few. A unique transmission function should be determined for all regularly used combinations of analyser settings to ensure consistency across multiple experiments. The transmission function will drift over time and when a significant change is made to the XPS instrument (e.g. after a bake-out), so it follows that without regular intensity calibration, the accuracy of analyses from a given spectrometer can become unreliable. This issue is compounded when comparing results from different laboratories – if these laboratories do not calibrate their analysers regularly, or are calibrating with respect to different

reference spectra, then comparing XPS data between laboratories will not be possible. In an era where the reproducibility of results in scientific publications is being called into question,⁶⁻⁸ it seems appropriate to highlight the calibration of the instruments that generate scientific data. Several successful international standards have already been published under the auspices of ISO/TC201/SC7 Electron Spectroscopy Standards. For example, there already exists an international standard for the calibration of energy scales in XPS (ISO 15472:2010),⁹ as well as the repeatability and constancy of the intensity scale (ISO 24237:2005).¹⁰ As of writing, there is no international standard protocol for calibrating the intensity scale of XPS spectrometers and this study is intended to both motivate and inform the generation of such standards. Many methods of intensity calibration and reference materials have been used for XPS instruments,¹¹⁻¹⁸ including the use of peak intensities,^{13,14,18} sample bias,¹⁴ elastic electron scattering from an isotropic electron gun,¹⁵ and general assumptions regarding the transmission of the analyser, to name a few.^{16,17} An intensity calibration method developed at NPL by Seah and Smith, compares experimental survey scans from copper, silver and gold to accurate reference spectra suitable for the instrument geometry.¹⁹⁻²¹ The ratio of the experimental result to the accurate reference spectra provides a representation of T with additional experimental noise, denoted Q . The rapidly changing intensities near the peaks are excluded for two reasons. Firstly, due to differences in instrument resolution, the widths and shapes cannot easily be adapted to provide a suitable match. Secondly, the large intensities at the peak are more likely to be in error due to detector saturation. An alternative method of determining a spectrometer's transmission function uses low-density polyethylene (LDPE, $(C_2H_4)_n$) as a secondary reference material (i.e. where the secondary reference spectrum is obtained from a spectrometer calibrated using a primary reference material such as gold, silver, and/or copper).²² LDPE is a thermoplastic widely used for containers, packaging, and tubing due to its excellent chemical resistance, flexibility, and toughness, as well as being non-toxic. If measured using XPS without preparation, LDPE will exhibit an O 1s peak on the surface due to long-term environmental oxidation and oxygenated functional groups in adventitious carbon contamination. Previous studies with Argon cluster ion sputtering show that this oxidation extends hundreds of nanometres into the material.²² It was confirmed that scraping an LDPE surface several times with a clean scalpel blade results in an oxygen-free surface which is insensitive to carbon contamination,²³

meaning that LDPE can be prepared *ex-situ* without the need for an ion source for sputtering. LDPE's chemical structure results in a spectrum that contains only carbon-related peaks (i.e. *C KLL* Auger feature, the *C 1s* core level, and valence states), with most of the measured photoelectron intensity comprising of inelastic background. The simplicity of the inelastic background means that it can be described by a relatively straightforward mathematical function. This mathematical description of the LDPE inelastic background is easily reproduced, continuous, and noise-free, which is a significant advantage over precious metal reference spectra using this calibration approach.

Under the auspices of Technical Working Area 2 (Surface Chemical Analysis) of the VAMAS organisation, an inter-laboratory study was conducted to further investigate XPS intensity calibration using LDPE as a reference material. The general aim of the study was to assess the validity of this calibration method across a wide range of laboratories with different XPS instruments. To date, this intensity calibration method has only been applied to XPS spectrometers which use a quartz-monochromated *Al K α* X-ray source. We did not receive any datasets from participants with non-monochromated X-ray sources or other anode targets such as Mg or Cr.

II. VAMAS PROJECT A27

A. Experimental

A 150 mm \times 150 mm \times 1 mm sheet LDPE (CAS No. 9002-88-4) was purchased from Goodfellow Cambridge Ltd (Huntingdon, UK) and cut into 50 individual pieces of approximately 10 mm \times 20 mm in size. The edges of the LDPE were trimmed with a clean scalpel to remove any ragged edges. Each piece of LDPE was then wrapped in aluminium foil to prevent degradation due to light until the sample was ready to be analysed by the participant. Two sheets of gold (AU153709, 99.99 % purity, 50 mm \times 50 mm \times 0.125 mm) were purchased from Advent Research Materials (Oxford, UK) and cut into 50 individual pieces of approximately 10 mm \times 10 mm in size. A plastic protective coating was attached to one side of each gold piece, which was removed using tweezers. Each piece of gold was cleaned by sonication in acetone for 15 minutes, then sonication in isopropanol for 15 minutes, followed

by a blow-dry with argon gas. Clean and sterile disposable scalpels (10A, Ref 0502) were procured from Swann-Morton Ltd (Sheffield, UK). The samples of LDPE and gold were placed into clean 1" Fluoroware wafer shipping containers (Megatech Limited, Huntington, UK). The samples were shipped to the participants along with a disposable scalpel and a hard copy of the protocol (see Supporting Information S1).

Participants were asked to mount the LDPE and gold samples onto their instrument's sample holder using any appropriate method. To prepare a fresh oxygen-free surface on the LDPE, participants were required to use the disposable scalpel to scrape the LDPE surface several times until it turned from shiny to matte in appearance. Participants were instructed to ensure that the direction of the last few scrapes were parallel to the X-ray source-analyser plane, to reduce X-ray shadowing effects from the topography generated on the LDPE during preparation. Participants were advised to prepare the LDPE surface immediately before introduction to UHV. If the procedure was performed correctly, then the O 1s signal in XPS should be below the detection limit which is less than 0.03 at%.²³ Similarly, the gold surface must also be contaminant-free, so participants were recommended that Ar-ion sputtering *in situ* be used. If they did not have an argon-ion sputtering gun, then they were asked to clean the gold surface using any method available in their laboratory, or to evaporate gold in situ. Again, all details on handling of samples can be found in the full protocol in the supporting information (See S1) or online.²⁴

B. Analysis

Prior to any analysis on the samples, participants were instructed to obtain an XPS survey spectrum without a sample in the analysis position and with the X-ray source switched off in order to obtain the 'dark noise' count of the XPS instrument. Participants were then asked to switch on their X-ray source and allow at least 30 minutes for the instrument to equilibrate. On the LDPE, participants were directed to select an area of analysis that was free of oxygen or any other contaminating species. If they were unable to find a clean area, then the LDPE should be scraped with the disposable scalpel again to generate a clean surface. LDPE is naturally insulating so participants were asked to use a low-energy electron source to compensate for surface charging, such that a single-component C 1s peak was

observed between 1200 eV and 1206 eV kinetic energy (for Al K α X-ray sources). Once the sample height and charge compensation settings were optimised, participants were then instructed to acquire three spectra on the LDPE: an initial full range survey, I_{initial} (180 eV to 1500 eV); a high kinetic energy (HKE) region spectrum, I_{HKE} (1195 eV to 1500 eV); and a final full range survey, I_{final} (180 eV to 1500 eV). The inelastic background in the HKE region of a LDPE spectrum has a low intensity compared to the rest of the spectrum. Participants were therefore asked to increase the acquisition time of I_{HKE} by a factor of ~30 to 40 by increasing the dwell time or by taking multiple acquisitions. This procedure was to be repeated on a fresh area of LDPE for each lens mode, pass energy, and slit combination. Participants were free to submit as many datasets for different analyser parameters as they wished. The minimum requirement was to submit at least one set of spectra acquired using a high pass energy (e.g. ≥ 90 eV) and a fully open entrance slit. For every set of analyser parameters used to acquire LDPE spectra, participants were also instructed to acquire a full survey spectrum (180 eV to 1500 eV) of sputter-cleaned gold using the same acquisition parameters as the full LDPE survey spectra. The full procedure for acquiring XPS spectra for the inter-laboratory study is described in the protocol in the supporting information (see S1). A copy can also be obtained from the NPL website.²⁴

All spectra were to be exported without any existing transmission function correction and the photoelectron intensity had to be reported in counts per second. Participants were provided with an MS Excel file with which to report their data and metadata back to NPL. This Excel file also provided a step-by-step procedure of how to generate a transmission function from LDPE, and participants were encouraged to attempt the analysis detailed in S1 themselves, although this was not a necessity.

C. Participants

Sample packages were distributed to 44 participating laboratories: 17 academic institutions, 16 industry partners (including 5 XPS instrument manufacturers), and 11 national facilities. At the conclusion of the study, 35 participants returned at least one dataset to NPL. Some participants provided multiple datasets for different instruments and operating conditions. Table I shows a list of participant codes with the corresponding details and geometry of their XPS instruments. A wide variety of instruments

from different manufacturers with different geometries allowed the applicability of the LDPE intensity calibration technique to be investigated thoroughly.

TABLE I. A list of participants' instruments with corresponding geometry. The terms a , b , and G are detailed in the results section and reference 25. Where angle b is not reported (—), an angle of 180° has been used to calculate G .

Code	Spectrometer Manufacturer / Model	Angle a	Angle b	G (2 s.f.)
AA	Thermo Scientific / Alpha 110	41	180	-0.52
AB	Kratos Analytical / Nova	54.7	180	-0.22
AC	PHI / VersaProbe II	45	180	-0.43
AD	Thermo Scientific / Theta 300	56	45	-0.11
AE	PHI / VersaProbe II	45	180	-0.43
AF	Kratos Analytical / Axis Ultra DLD	60	180	-0.10
AG	Thermo Scientific / Nexsa	60	180	-0.10
AH	PHI / VersaProbe II	45	180	-0.43
AI	Scienta Omicron / UHV System	54.7	—	—
AJ	ULVAC-PHI / Quantum 2000	45	180	-0.43
AK	SPECS / PHOIBOS 150	54	180	-0.23
AL	PHI / VersaProbe II	45	180	-0.43
AM	Kratos Analytical / Axis 165	54	180	-0.23
AN	SPECS / PHOIBOS 150	55	180	-0.21
AO	Kratos Analytical / Axis Ultra DLD	60	180	-0.10
AP	Kratos Analytical / Axis Ultra DLD	60	180	-0.10
AQ	Thermo Scientific / VG ESCALAB 250 Xi	58	180	-0.14
AR	PHI / VersaProbe III	45	180	-0.43
AS	Kratos Analytical / Axis Ultra DLD	60	180	-0.10
AT	SPECS / PHOIBOS 150	55	180	-0.21
AU	PHI / Quantera II Hybrid	45	180	-0.43
AV	Surf Sci / SSX-100	75.5	120	0.37
AW	Kratos Analytical / Axis Supra DLD	60	180	-0.10
AX	SPECS / PHOIBOS 150	55	180	-0.21
AY	ULVAC-PHI / Quantera SXM	45	180	-0.43
AZ	Thermo Scientific / Nexsa	60	180	-0.10
AZ	Thermo Scientific / VG ESCALab 250Xi	58	180	-0.14
BA	Kratos Analytical / Axis Supra	60	180	-0.10
BB	PHI / Quantera Hybrid	45	—	—
BC	ULVAC-PHI / Quantera SXM	45	180	-0.43
BD	SSI / S-Probe	55	—	—
BE	Kratos Analytical / Axis Ultra DLD	60	180	-0.10
BF	Thermo Scientific / ESCALAB Xi+	58	180	-0.14
BG	Thermo Scientific / VG ESCALAB 220i-XL	54.7	180	-0.22
BG	ULVAC-PHI / ESCA 5800	60	54.7	0.01
BH	Kratos Analytical / Axis Ultra DLD	60	180	-0.10
BI	Kratos Analytical / Axis Ultra DLD	60	180	-0.10

III. RESULTS AND DISCUSSION

For each participant and analyser settings, the set of LDPE spectra, I_{initial} , I_{HKE} , and I_{final} were combined into a single spectrum in counts per second and then a dark noise correction, D , was subtracted. The combined spectrum for a given participant is denoted I_{PE} herein. Equation (4) in the “XPS Intensity Calibration Protocol using LDPE” section describes how I_{PE} is calculated. The corresponding sputter-cleaned gold spectrum, acquired using the same acquisition parameters, is denoted I_{Au} . Each participant's I_{PE} (and I_{Au}) was visually checked for any sign of contamination and for evidence of internal scattering in the analyser; any factors in the participants' data which may affect the final transmission function were noted in Table S2 in the supporting information. By dividing I_{PE} by the LDPE reference spectrum, a noisy representation of the XPS instrument's intensity response, Q , is obtained, from which the instrument's transmission function, T , can be determined. However, this is not the full story.

A. Geometry-corrected Reference Spectra for Polyethylene

In the original publication²² and VAMAS protocol²⁴ (S1), the reference LDPE spectrum used for intensity calibration is generated using a six-component mathematical description which accounts for all of the observed features in the LDPE background. The fitting parameters used to reproduce the LDPE reference spectrum were determined by fitting an NPL transmission function-corrected LDPE spectrum acquired using the NPL spectrometer (Kratos Axis Ultra DLD). The resulting fitting parameters and background description accurately described the measured LDPE spectrum with a systematic error within 1 % across most of the spectrum, increasing to 5 % around the region of the *C KLL* features. However, these fitting parameters are only valid for instruments with a similar geometry to the NPL spectrometer where the angle between the incoming X-ray vector and the sample-to-analyser vector, a , is 60° , and the angle between the anode-monochromator-sample plane and analyser-sample-monochromator plane, b , is 180° . For this work, the geometry of a given instrument is termed $[a, b]$, so the NPL spectrometer has a geometry of $[a, b] = [60^\circ, 180^\circ]$. As evident in Table I, a variety of instrument geometries were provided, so it was necessary to investigate the effects of

instrument geometry on the LDPE spectrum. This is because Al $K\alpha$ X-rays that interact with a quartz crystal monochromator have some degree of partial polarisation which alters the distribution of photoelectrons as a function of the instrument geometry.²⁶ A new mathematical description of the LDPE reference spectrum for all geometries has been developed and disseminated by Shard & Reed.²⁵ It uses a practical form of the angular distribution of photoelectron emission due to partially polarised X-rays and angular anisotropy in photoemission, G , and the energy-dependent geometry correction for elastic scattering, f . Herein, the LDPE reference spectrum for a specific instrument geometry is denoted as I_{ref} . As shown in Table I, the most common geometries encountered during the VAMAS study were $[45^\circ, 180^\circ]$ and $[60^\circ, 180^\circ]$ for which $G = -0.43$ and -0.10 , respectively. Suitable contaminant-free I_{PE} and I_{Au} were selected from the participants that had instrument geometries of $[45^\circ, 180^\circ]$ and $[60^\circ, 180^\circ]$. The I_{PE} spectra were intensity calibrated using a transmission function determined using their corresponding I_{Au} and the NPL reference spectra for gold. Then, an average intensity-calibrated I_{PE} with intensity units I_X (consistent with the NPL calibration method) was determined for each geometry. Using these averaged spectra, shown in Figure 1A, and the methodology described by Shard & Reed,²⁵ I_{ref} can be determined for any $[a, b]$ instrument geometry. Figure 1B shows I_{ref} for varying angles of a and b . Figure 1C demonstrates why this geometry correction is required: when I_{PE} from participant AA is divided by the original LDPE reference spectrum described by Shard & Spencer,²² a significant intensity difference is observed between the high and low kinetic energy sides of the C 1s region. This is because the original LDPE reference spectrum is only valid for an instrument geometry of $[60^\circ, 180^\circ]$, and the background intensity dependence either side of the C 1s region does not match the I_{PE} from participant AA. Using an I_{ref} which is geometry corrected for $[41^\circ, 180^\circ]$, the resulting ratio $Q = I_{\text{PE}}/I_{\text{ref}}$ does not exhibit the intensity difference around the C 1s region and, therefore, a more sensible description of the spectrometer's transmission function is obtained.



This is the author's peer reviewed, accepted manuscript. However, the online version of record will be different from this version once it has been copyedited and typeset.
PLEASE CITE THIS ARTICLE AS DOI: 10.1116/6.0000577

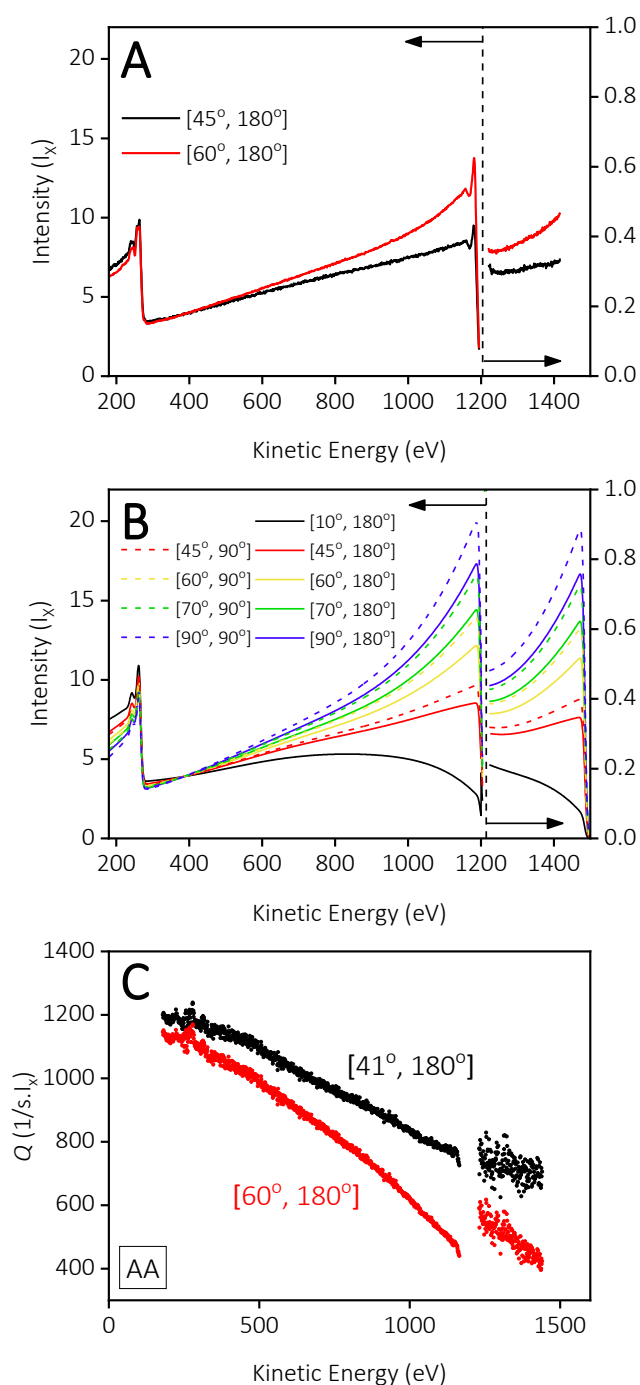


FIG. 1. Geometry correction to the LDPE I_{ref} . In panels (A) and (B), the black dashed line denotes an axis break in kinetic energy where the low KE region scale is on the left, and the high KE region scale is on the right. (A) The two calibrated LDPE spectra from the VAMAS participants data. The [45°, 180°] spectrum is shown in black, and the [60°, 180°] spectrum is shown in red. (B) Geometry-corrected LDPE I_{ref} calculated using the method described by Shard & Reed.²⁵ A reference spectrum for any combination of a and b can be calculated. (C) The effect of using the new geometry correction on Q_{PE} for participant AA. The red dots show the values of Q_{PE} when calculated using a LDPE I_{ref} for the [60°, 180°] geometry, and the black dots show the values of Q_{PE} when the correct [41°, 180°] geometry LDPE I_{ref} is used.

Figure 2A shows examples of I_{PE} from several participants. For the purposes of data visualisation, these spectra have been normalised to the intensity of the $C\ KLL$ feature, I_{C-KLL} , in order to compare the inelastic background intensities on the same scale. The corresponding sputter-cleaned gold spectra, I_{Au} , (also normalised to the $Au\ NNV$ area, I_{Au-NNV} , for comparison) acquired using the same analyser and acquisition parameters are shown in Figure 2B. The significant variability of background shapes and intensities of I_{PE} and I_{Au} between different participants, visually demonstrates how variable the responses of different instruments are, and therefore how important the transmission function correction is. I_{PE} was then divided by the corresponding I_{ref} to provide a representation of the transmission function Q . Excluding the regions around the $C\ 1s$ peak and valence band, the values of Q can then be fitted using the same functional form used in the protocol (see equation 3 in Shard & Spencer²²) in order to determine a noise-free description of the transmission function T . The values of Q and the functions representing T , shown in Figure 2C and Figure 2D, are calculated using I_{PE} and I_{Au} shown in Figure 2A and Figure 2B respectively. Herein, for each dataset from every participant, the values of Q determined from I_{PE} and I_{Au} , are denoted as Q_{PE} and Q_{Au} respectively. Similarly, the functions of T determined from I_{PE} and I_{Au} , are T_{PE} and T_{Au} respectively. The discrepancies between T_{PE} and T_{Au} for a given participant are discussed later in this report. Once T has been determined for a given set of analyser parameters, then any spectrum acquired using the same parameters can be intensity-corrected using the same T . Figure 2E and Figure 2F show the initial I_{PE} and I_{Au} divided by their respective T . As expected in the case of the LDPE spectra, the intensity corrected spectra take on the shape and the intensity of I_{ref} . Participants AX, AT, and AI acquired I_{PE} using the same $[55^\circ, 180^\circ]$ geometry but the shape and intensity of the spectra across the kinetic energy range vary significantly. After intensity-correction using their corresponding instrument's T_{PE} , all the LDPE spectra acquired using the $[55^\circ, 180^\circ]$ geometry are identical. Similarly, after intensity-correction using T_{Au} , all the gold spectra exhibit nearly identical shapes and intensities; the relative effect of instrument geometry is less pronounced due to the larger average intensity of gold compared to LDPE. Figure 2E and Figure 2F demonstrate that the LDPE method of calibration can reach the same result as the NPL intensity calibration method using gold (plus silver and copper) reference spectra. Quantitative analysis of the same sample across these XPS instruments would provide consistent results in terms of peak intensities now that they have functions



This is the author's peer reviewed, accepted manuscript. However, the online version of record will be different from this version once it has been copyedited and typeset.
PLEASE CITE THIS ARTICLE AS DOI: 10.1116/6.0000577

of T that are traceable to the same set of reference spectra. The values of T are valid for as long as the analyser is stable (i.e. the transmission function may drift over time as the analyser ages) and no major changes take place in the XPS instrument (e.g. venting, bake out, maintenance on components of the XPS instrument, etc.).

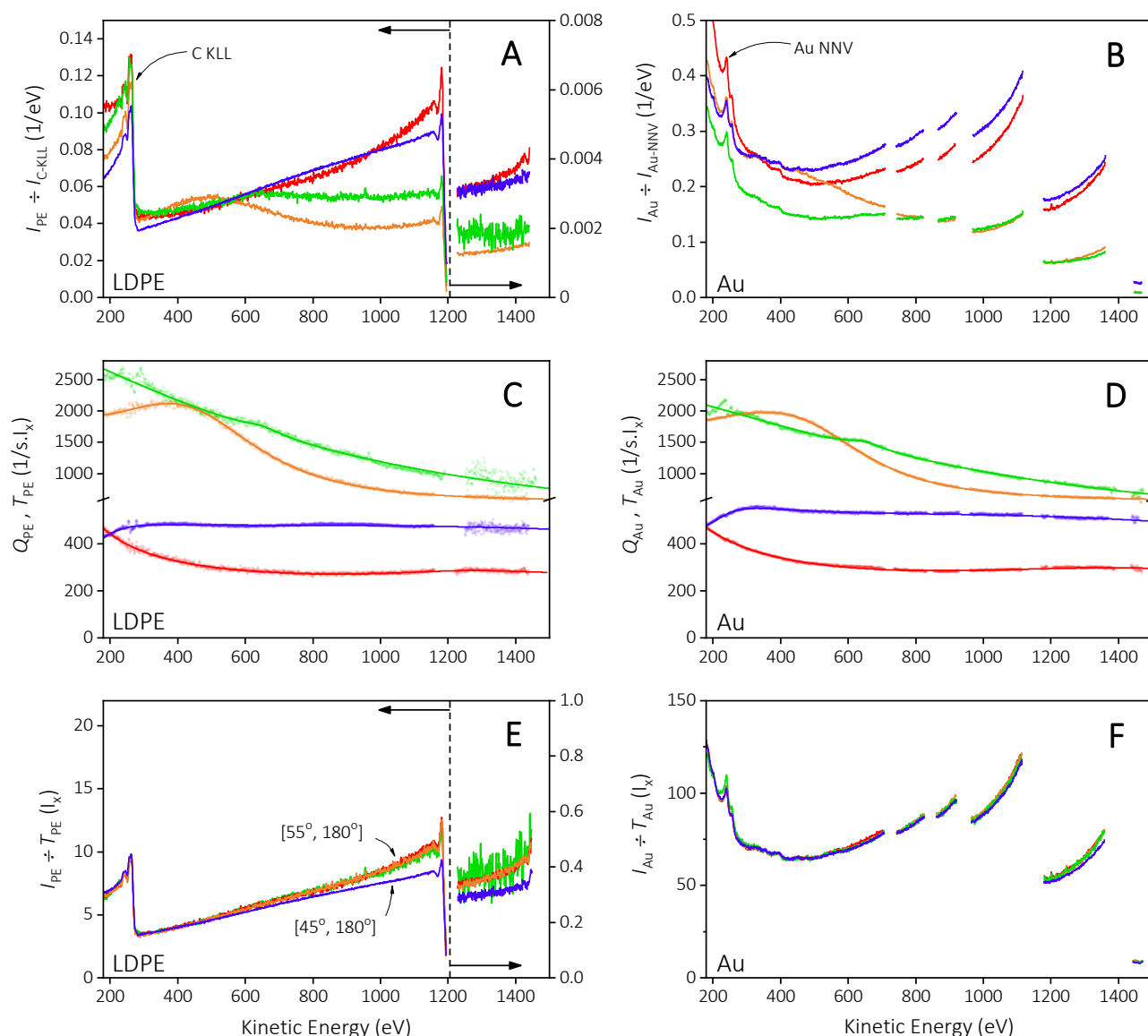


FIG. 2. Transmission function correction using LDPE and Au spectra from participants AX (red), AJ (blue), AI (green), and AT (orange). **(A, B)** LDPE and Au spectra normalised to the C KLL and Au NNV Auger peaks respectively. **(C, D)** Q (dots) and T (lines) calculated using LDPE and Au intensity calibration methods respectively. **(E, F)** Intensity-corrected LDPE and Au spectra using T_{PE} and T_{Au} respectively.

B. Analysis of VAMAS Data

As discussed previously,²² the use of the original LDPE reference spectrum for intensity calibration resulted in transmission functions with absolute values within 10 % of those calculated using the NPL software and systematic deviations of less than 5 %. However, error analysis was only performed on one dataset from the NPL spectrometer with a [60°, 180°] geometry. In line with the aims of VAMAS study A27, the applicability of the LDPE intensity calibration method depends on how this error varies from instrument to instrument and when different analyser settings are used. T_{Au} is directly traceable to the NPL gold reference spectrum and so, assuming that the gold foil was sputter cleaned correctly and the spectra did not exhibit any scattering or non-linearity effects, T_{Au} can be used to assess any absolute offset and systematic deviations in T_{PE} . There are limitations to this assessment, which will be discussed later, but it can at least show whether the LDPE intensity calibration method reaches a similar transmission function as gold. The error ratio, R_i , is equal to T_{PE} divided by T_{Au} at kinetic energy data point i . If T_{PE} has the exact same intensity of T_{Au} , then R_i will be equal to unity across all kinetic energies. In practice, R_i deviates from unity across all energies in all the participants' datasets, so it helps to define terms which characterise these deviations. The average offset factor, Δ , of T_{PE} with respect to T_{Au} can be determined by subtracting unity from the mean value of the error ratio, \bar{R} , over all values of i up to the number of kinetic energy data points, n , that is:

$$\Delta = \left[\left(\frac{1}{n} \sum_{i=1}^n R_i \right) - 1 \right] = [\bar{R} - 1] \quad (2)$$

The systematic deviation, Σ , which characterises the difference in the shape of T_{PE} with respect to T_{Au} is calculated as twice the standard deviation of R_i across all values of i , i.e. $2\sigma_R$, divided by \bar{R} , such that:

$$\Sigma = \frac{2\sigma_R}{\bar{R}} \quad (3)$$

The confidence interval of Σ , based on an interval of $\pm 2\sigma$, is 95%, which is enough given that the accuracy of quantification for photoelectron peaks is further limited by errors in the inelastic background selection and effective attenuation length. For clarity, a graphical representation of R_i , Δ ,

and Σ from participant AM is shown in the supporting information (Figure S7). In this work, we calculated the average and standard deviation of R_i over a range of 180 eV to 1500 eV kinetic energy. The values of Δ and Σ are relative uncertainties and so percentage values are obtained by multiplying them by 100. The percentage values of Δ (%) and Σ (%) for each dataset provided by the participants are plotted against each other in Figure 3. The black dashed lines in Figure 3 show the mean value of Δ (%) and Σ (%) when considering all the participants' data, which correspond to values of 3.0 % and ± 6.5 % respectively (1 d.p.). Comparing the datasets with respect to parameters such as lens modes, slit width, and certain other instrumental factors is not possible as the design and operation of each participant's analyser is not necessarily comparable. However, the analyser pass energy is a measurable numerical parameter and has a predictable effect on the participants' data, i.e. generally the larger the pass energy, the larger the intensity measured on the detector. The reported pass energies from the participants are separated into one of two groups: low (< 90 eV), and high (≥ 90 eV). In Figure 3, the blue dashed lines represent the mean value of Δ (%) and Σ (%) for datasets acquired at low pass energies, and the red dashed lines represent the mean value of Δ (%) and Σ (%) for datasets acquired at high pass energies. For low pass energies, the mean values of Δ (%) and Σ (%) are -4.9 % and ± 8.8 % respectively, and for high pass energies, the mean values are 5.8 % and ± 5.7 %. The average results of Δ (%) and Σ (%) for low and high pass energies are reported in Table II.

TABLE II. Average values of Δ (%) and Σ (%) for all participants' datasets, low pass energy datasets, and high pass energy datasets.

Dataset Grouping	Δ (%)	Σ (%)
All participants	3.0 %	± 6.5 %
Low pass energy (< 90 eV)	-4.9 %	± 8.8 %
High pass energy (≥ 90 eV)	5.8 %	± 5.7 %

The average offset remains within 10 % for all the participants' datasets, as was similarly reported in the original NPL publication on LDPE calibration. However, from the point of view of obtaining correct relative intensities, the systematic deviation parameter is more important, and the effect of pass energies on Σ (%) is more pronounced. We observe that the lower pass energy group has more systematic

deviation compared to the high pass energy. Increasing the analyser pass energy increases the intensity of the photoelectron signal whilst sacrificing energy resolution, so in general we observe greater signal-to-noise across all kinetic energies in the high pass energy scans. The smaller signal-to-noise in the low count rate, low pass energy scans seems to result in more noise in the Q_{Au} and Q_{PE} values which inevitably results in more variability in the fitting between T_{Au} and T_{PE} .

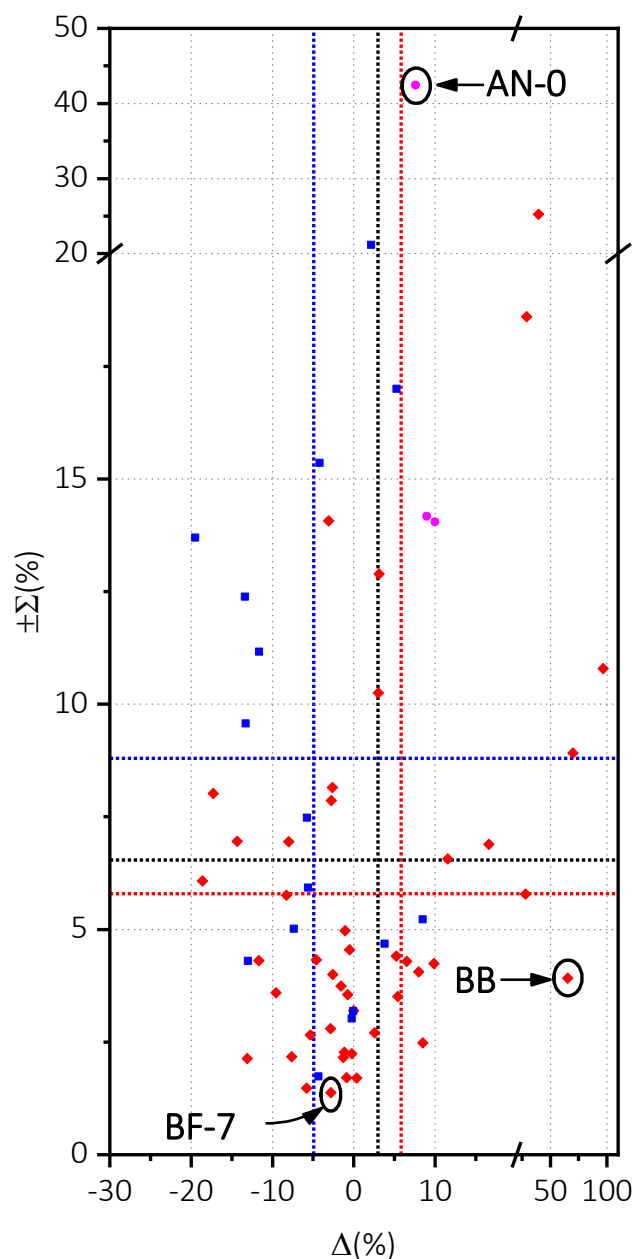


FIG. 3. The values of $\Delta(\%)$ and $\pm\Sigma(\%)$ for all the participants' datasets, where the red diamonds (\blacklozenge) refer to high pass energy datasets (≥ 90 eV), and the blue squares (\blacksquare) refer to low pass energy data sets (< 90 eV). The dashed lines show the average values of $\Delta(\%)$ and $\pm\Sigma(\%)$ from all datasets (black), high pass energies (red),

and low pass energies (blue). The pink dots (●) refer to datasets which have been replaced with newer data from the same participant, and do not contribute to the average values of Δ (%) and $\pm \Sigma$ (%).

Three datapoints in Figure 3 have been highlighted which exhibit either a large Σ (%) (AN-0), a large Δ (%) (BB), or both a low Σ (%) and Δ (%) (BF-7). The values of Q and T derived from LDPE and gold for these three datapoints are shown in Figure 4 for demonstrative purposes. Similar figures have been produced for each set of the data submitted by the participants and can be viewed in the supporting information (Section S2). Participant AN-0 has produced a large Σ (%) due to large differences in the slope and shape between their T_{Au} and T_{PE} . Assuming a lack of operator error or sample issues such as contamination, this case is a cause for concern as the disparity between the LDPE and gold data could suggest an underlying issue with the spectrometer. Participant BB has produced a reasonable Σ (%) (in fact, better than the average Σ (%) for high pass energies), but there is a significant offset between T_{Au} and T_{PE} resulting in a large Δ (%). Although not ideal, this situation is workable as it still permits transmission function correction with an aim to extract relative peak intensities for quantification. The reason for the offset is also easier to diagnose; the most common cause for the offset between T_{Au} and T_{PE} is either a difference in the X-ray source power, an insufficient warm up time for the X-ray source between the LDPE and gold scans, or failure to report data in counts per second. There is also supporting evidence to suggest that the topography of the LDPE after scraping with a scalpel can influence the absolute intensity of the LDPE spectrum (Figure S8); differences in Δ (%) up to 20 % are possible through excessive scraping, but the sample orientation after scraping is not a critical factor. The small values of Δ (%) and Σ (%) obtained from participant BF-7 present a case where T_{Au} and T_{PE} are nearly identical, showing that, in the absence of instrumentation issues and operator error, the LDPE and gold intensity calibration methods can produce the same result.

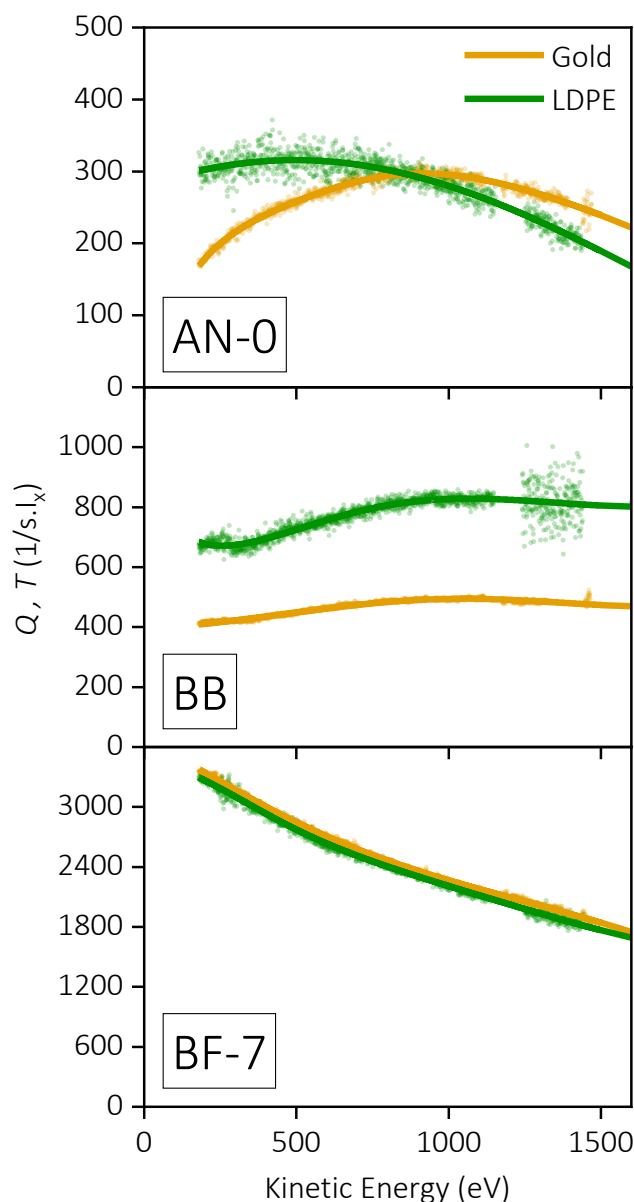


FIG. 4. The Q (dots) and T (lines) calculated from LDPE (green) and Au (orange) respectively for participants AN-0, BB, and BF-7, which have been highlighted in Figure 3.

In the original protocol, participants were instructed to acquire a set of ‘dark noise’ spectra with the X-ray source on and off. The average counts per second of this dark noise, D , would be subtracted from I_{PE} . This step is important, especially for the low count rate region on the low kinetic energy side of the C 1s peak. If the dark noise intensity is not removed from I_{PE} , then a significant intensity step is observed in Q_{PE} across the C 1s region, which cannot be fitted by the functional form of T_{PE} . In most cases, an instrument’s dark noise will be constant across the entire kinetic energy range, and so dark noise correction is relatively straight forward. However, in the study we encountered instances where the dark noise is dependent on the kinetic energy, and needed to be given a functional form, $D(E)$. Such a case

is not ideal as it means that there is an underlying issue with the spectrometer. Figure 5A shows an example, from this VAMAS study (participant AX), of a spectrometer exhibiting a variable dark noise intensity. Assuming the necessary troubleshooting has been attempted and $D(E)$ is stable, then some arbitrary, noiseless function $f(E)$ such as the one used in Figure 5A can be used to fit $D(E)$ which can then be used to perform the dark noise correction. It should be noted that the function $f(E)$ is not a valid physical model describing the response of the spectrometer, and only serves to correct the dark noise intensity. The effect of the variable dark noise correction on Q_{PE} (and, hence, T_{PE}) is shown in Figure 5B, where the LDPE intensity calibration method has been applied to I_{PE} with *i*) constant D correction (black), and *ii*) variable $D(E)$ correction (red). In both cases, the values of Q_{PE} are similar on the low kinetic energy side of the C 1s region. However, in the case of constant dark noise correction, there is an intensity step in Q_{PE} over the C 1s region which is caused by the mismatch of true dark noise intensity between the low and high kinetic energy regions. By using $f(E)$ to correct the variable dark noise, the intensity mismatch around the C 1s region in Q_{PE} is drastically reduced, and the fitted T_{PE} more pertinently describes the shape of Q_{PE} between 1000 eV and 1500 eV, which accords more closely with T_{Au} . Whilst the agreement between the resulting T_{PE} functions using different dark noise corrections is $< 1\%$ at lower kinetic energies, there is a maximum difference of 4.9 % at higher kinetic energies. This, therefore, demonstrates the importance of the correct subtraction of dark noise intensity from I_{PE} . Whilst the variable dark noise correction is valid in this case, it should be stressed that such an issue with the spectrometer should be investigated prior to the calibration procedure.

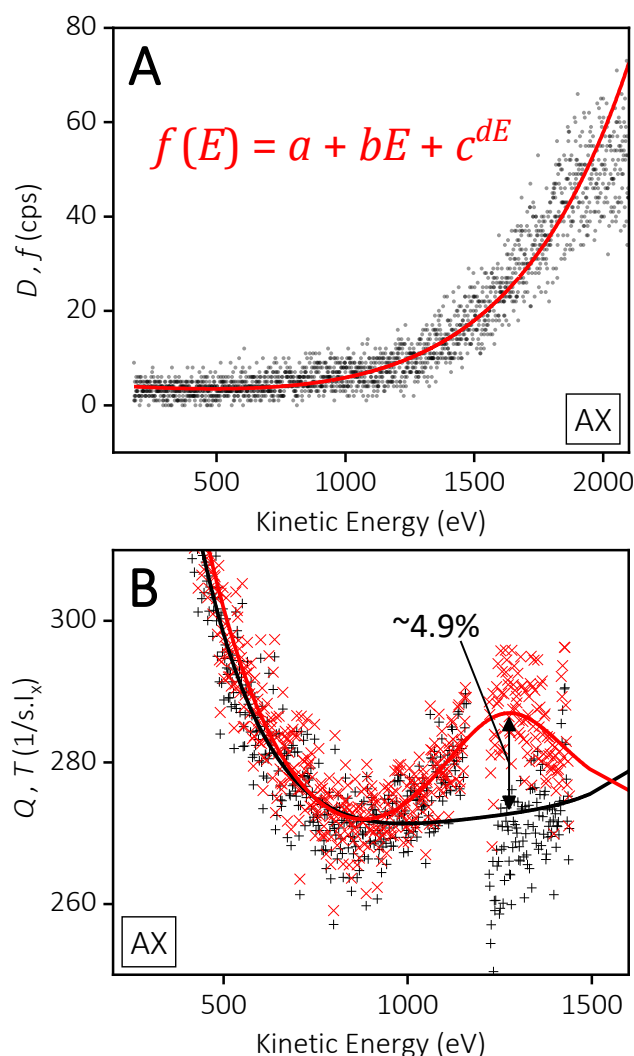


FIG. 5. (A) The variable dark noise $D(E)$ recorded by participant AX (black dots) fitted with an arbitrary function $f(E)$ (red line). (B) The effect of correcting $D(E)$ shown by participant AX. The difference in Q_{PE} when corrected with a constant dark noise value (black +) or $f(E)$ (red x). In the high kinetic energy, low count rate region above the C 1s, the difference in Q_{PE} is $\sim 4.9\%$.

Some participants' instruments exhibited behaviour that could not be rectified within the scope of this inter-laboratory study. In these cases, underlying instrument issues caused large deviations in Q which could not be captured by the functional form for T . Figure 6 shows where unknown issues in participant BD's analyser caused large deviations in Q_{PE} resulting in an intensity mismatch across the C 1s region. Furthermore, a significant loss of counts occurs in the low kinetic energy region (i.e. < 300 eV) for both Q_{PE} and Q_{Au} . These features meant that the functional form of T could not converge on Q_{PE} and so the fitting failed. Whilst not an ideal outcome, the low count rate of the LDPE spectrum clearly highlights an underlying issue with the spectrometer, which the gold spectrum did not pick up on. Hence the LDPE

intensity calibration method has potential as a tool to diagnose subtle issues in the analyser which manifest at low count rates.

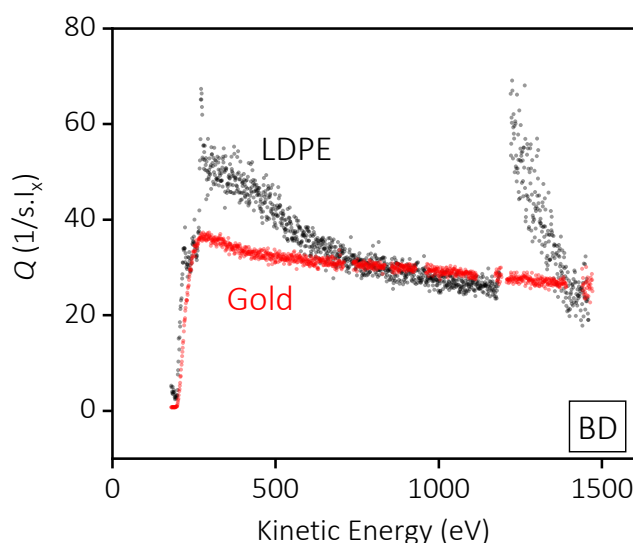


FIG. 6. Q_{PE} and Q_{Au} from participant BD showing significant internal scattering and intensity loss in the low kinetic energy region.

As previously mentioned, LDPE can be cleaned *ex-situ* using nothing more than a clean scalpel or razor. It is important that the LDPE surface is scraped sufficiently to remove any surface oxidation or contamination, as obtaining meaningful values of Δ (%) and Σ (%) is highly dependent on acquiring both the LDPE and gold spectra under optimum conditions. Figure 7 demonstrates the effect of surface contamination on Q_{PE} and T_{PE} . The I_{PE} from participant AD-2 contained a significant O 1s peak and O KLL Auger feature, as well as smaller silicon and sulphur peaks (2p and 2s) in the HKE region. The inelastic background of the contaminant peaks is not accounted for by the LDPE reference spectrum and, thus, the resulting Q_{PE} contains peaks and inelastic background steps related to these elements. The increase in the inelastic background due to the silicon and sulphur peaks result in an intensity mismatch across the C 1s region which T_{PE} cannot account for. The cumulative contributions of the background from the contamination peaks also result in a large deviation in Q_{PE} compared to Q_{Au} at low kinetic energies. Usually this problem can be easily resolved by removing the LDPE sample from the XPS instrument and repeating the sample preparation step in the procedure. However, persistent contamination that cannot be removed may suggest either a bulk contaminated LDPE sample (requiring the LDPE to be replaced) or a contaminated XPS instrument (requiring a bake out). Another source of

contamination to consider is the scalpel or razor used in the LDPE preparation, for example commercial razor blades often have a polytetrafluoroethylene coating which could contaminate the LDPE surface with fluorine. In the case of gold, if there is any carbon contamination, this will cause a decrease in the count rate at lower kinetic energies due to the attenuation of the gold Auger features. Consequently, a significant deviation in Q_{Au} (and by extension, T_{Au}) at low kinetic energies, will result in a similar low kinetic energy deviation in R_i when compared to data from a clean oxygen-free LDPE spectrum. It follows then that the values of $\Delta(\%)$ and $\Sigma(\%)$ will not reflect a fair comparison, at least from a procedural perspective. Participant datasets with significant contamination on either the LDPE or gold have been identified in Table S2.

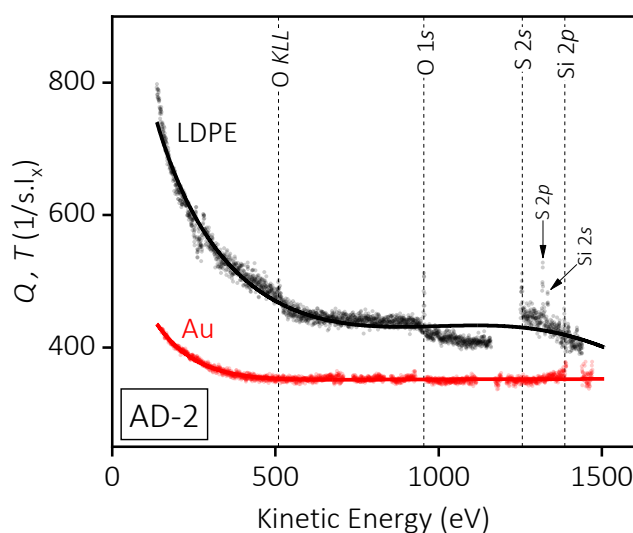


FIG. 7. Q_{PE} and Q_{Au} from participant AD-2 demonstrating the significant effect that contamination on the LDPE produces. Other elements (i.e. O, Si, and S) on the LDPE cause steps, discontinuities, and deviations in Q_{PE} .

In the VAMAS protocol, participants were instructed to obtain LDPE and gold spectra using the same analyser settings and X-ray source power. However due to the orders-of-magnitude difference in photoelectron yield between these samples, acquisition settings that produce an appropriate signal-to-noise for the LDPE may result in detector saturation when measuring the gold.²⁷ Due to the order of the original protocol, it is probable that participants measured the LDPE first, and then used the same parameters with the gold measurement; thus, resulting in detector saturation. Some participants anticipated this issue and submitted gold spectra acquired using a lower X-ray anode emission current to avoid saturation. In this situation, a correction factor equal to the ratio of the differing X-ray anode

emission currents is applied to I_{Au} or I_{PE} to reduce the X-ray power induced offset between Q_{Au} or Q_{PE} . Some participants' datasets exhibited evidence of detector non-linearity as a result of detector saturation when measuring the gold. Figure 8A shows Q_{Au} and T_{Au} from participant AF. There is an uncharacteristic intensity increase around the low kinetic energy side of the gold peaks that T_{Au} does not fit correctly. The participant's gold spectrum recorded a maximum Au $4f_{7/2}$ count rate of ~ 1.95 Mcps and an average background count rate of ~ 340 kcps. In comparison, I_{PE} does not exhibit any evidence of detector saturation, and the maximum C $1s$ count rate was ~ 300 kcps. Figure 8B shows the measured intensity of I_{Au} plotted against the gold reference spectrum I_{Au-ref} (geometry corrected) multiplied by T_{PE} (given that T_{PE} should not be affected by any non-linearity arising from detector saturation). The red dots show the actual gold intensity recorded by the participant and the black line shows the expected gold intensity (which has a slope of ~ 1.12 to correct for the absolute offset between T_{Au} and T_{PE}). If the detector has a linear response, then we would expect the red dots and black line to agree, but we observe a clear detector non-linearity in participant AF's data above ~ 300 kcps. Due to the low photoelectron yield from LDPE, non-linearity effects are unlikely to occur during the calibration procedure, especially using standard acquisition parameters – regardless, the non-linear regime of the detector should be determined prior to calibration in order to be sure.

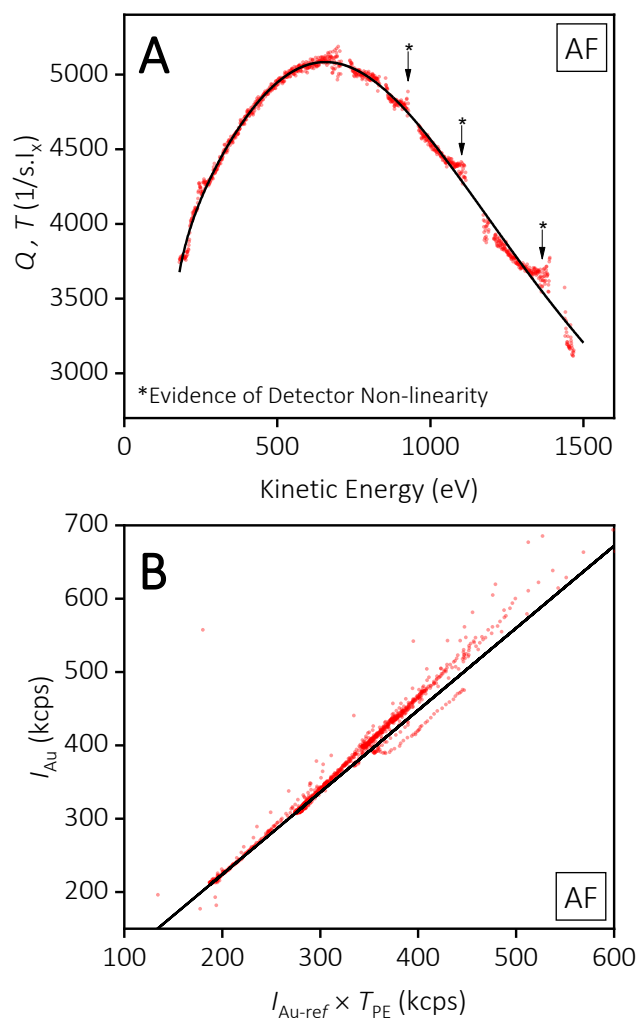


FIG. 8. **(A)** Q_{Au} (red dots) and T_{Au} (black line) from participant AF. The marked positions show an uncharacteristic increase in the spectrum intensity near the gold peaks which may affect the fitting of T_{Au} . **(B)** Plotting the actual I_{Au} from participant AF (red dots) against the expected I_{Au-ref} (intensity corrected using T_{PE} , and corrected for average offset) (black line) reveals detector non-linearity above 300 kcps.

This work has demonstrated the huge variability between the transmission functions of XPS instruments from around the surface analysis community, in different laboratories, and using different commercial and home-built system. At a minimum, it is important for XPS users to report the details of their instrument calibration, both in the energy and intensity scales, in their reports and publications in the interest of reproducibility.²⁸ However, we hope that this work has further communicated the vital need to use a consistent and traceable method of intensity calibration for the entire community of XPS users. This can be achieved via two avenues: *i)* Table I shows that most of the participants are using commercial XPS systems, so it is probably reasonable to say that most of the surface analysis (and surface science) laboratories in the community are operating commercial systems as well. As of writing,

XPS instrument manufacturers use their own in-house methods for intensity calibration, which may differ from company to company. An initial calibration is performed when the XPS instrument is sold, but this transmission function will quickly become redundant. If the user does not regularly calibrate their instrument, relying instead on the initial calibration file (assuming they are aware of it), then the results of their analyses will quickly become unreliable and irreproducible. An agreed traceable calibration method could be implemented by XPS instrument manufacturers and could either be integrated into both hardware and software allowing the user to easily and regularly calibrate the instrument, or be provided as part of regular engineer servicing and maintenance. This effort could be underpinned by *ii*) the publication of an international standard in collaboration with national metrology laboratories; As of writing, there is no international standard protocol for calibrating the intensity scale of XPS spectrometers. Such a standard should consider the alternative instrument configurations of its users, for example, some instruments do not have an *in-situ* sputtering source for cleaning sample surfaces. The LDPE intensity calibration technique can contribute to such an international standard, allowing XPS users that do not have a sputter source to calibrate their instruments. Furthermore, it is a method that does not require noisy experimental reference spectra (the raw data from which would need to be stored in an accessible database), instead using a well-defined mathematical description that is applicable to all instrument geometries.²⁵

Following the analysis of the inter-laboratory study data and feedback from the participants, an improved calibration protocol using LDPE has been detailed in the next section. This protocol addresses the issues brought up during the analysis with regards to geometry correction, variable dark noise correction, X-ray source warm up, sample contamination, analyser internal scattering, and detector non-linearity.

IV. XPS INTENSITY CALIBRATION PROTOCOL USING LDPE

The following intensity calibration protocol can be used for XPS instruments using a monochromatic Al K α X-ray source and a low-energy electron source to neutralise surface charging. This protocol is

applicable to instruments of any $[a, b]$ geometry. This protocol requires a sample of high-purity LDPE and a clean scalpel or razor blade; N.B.: avoid using commercial razor blades, as they are typically coated in a lubricating fluorocarbon which will contaminate the LDPE. It is advised that you read all the instructions carefully, ensuring that they are understood before continuing with the calibration procedure.

1. Determine the XPS instrument geometry, i.e. the source-to-analyser angle ' a ', and the dihedral angle between the monochromator-sample-analyser plane and the anode-monochromator-sample plane ' b '. See Shard & Reed.²⁵
2. Mount a piece of LDPE onto a sample holder. A single piece of LDPE may be used for many calibrations if it has been stored in a cool and dark place, away from sources of fouling (e.g. solvent fumes, heat, damp, prolonged light exposure, etc.). If the LDPE shows any sign of contamination or fouling, replace the LDPE stock with new material.
3. Using a clean scalpel or razor blade, firmly scrape the surface of the LDPE several times until the surface turns from shiny to matte in appearance. Do not over-scrape the LDPE; we recommend no more than 10 scrapes for a given area.
4. Immediately place the LDPE into the instrument load lock, and pump down. As soon as possible, transfer the LDPE sample into the analysis chamber.
5. With the LDPE sample positioned away from the analysis position, acquire a survey spectrum between 180 eV and 1500 eV kinetic energy to determine the 'X-rays OFF dark noise'.
6. Switch on the X-ray source and wait at least 30 minutes for the source to equilibrate. Some older or aging instruments may require longer.
7. With the LDPE sample positioned away from the analysis position, acquire a survey spectrum between 180 eV and 1500 eV kinetic energy to determine the 'X-rays ON dark noise'. Compare the dark noise surveys and ensure that they have a constant intensity across the entire spectrum range. If the dark noise varies across the survey range, this may be indicative of an instrument issue which should be investigated. As a temporary approach, prior to the instrument being

fixed, if the dark noise intensity follows a predictable and stable pattern, a suitable function $D(E)$ may be fitted to the dark noise intensity which can then be used.

8. Move the LDPE into the analysis position. Optimise the sample-to-analyser distance and the charge compensation settings to maximise the C 1s intensity. Ensure that the C 1s is a single peak positioned between 1201 eV and 1206 eV kinetic energy.
9. Acquire a survey spectrum between 180 eV and 1500 eV kinetic energy using a high pass energy (i.e. ≥ 90 eV). Ensure that no elements other than carbon are detected. If other elements are detected, move around the LDPE surface to find a contaminant-free area. If unsuccessful, remove the LDPE from the XPS instrument and start the calibration protocol again from step 3. If this fails again, check the LDPE, scalpel, laboratory, and XPS instrument for sources of surface contamination. If necessary, use a fresh piece of LDPE.
10. Change the spectrometer settings to the mode that requires calibration (i.e. pass energy, lens mode, exit and entrance slit, etc.) and allow enough time for the instrument to equilibrate if required. From this point forward, use the same acquisition parameters for all scans unless stated otherwise. Do NOT switch the X-ray source off between scans.
11. Acquire a survey spectrum between 180 eV and 1500 eV kinetic energy, E . Designate this spectrum $I_{\text{initial}}(E)$.
12. Acquire a survey spectrum between 1195 eV and 1500 eV kinetic energy with an acquisition time ~ 30 to 40 times longer than $I_{\text{initial}}(E)$. This can be done by either increasing the dwell time per step or increasing the number of acquisitions / sweeps / repetitions. Take an average of all spectra acquired in the HKE region in counts per second. (N.B. averaging XPS spectra needs to be done in counts per second, your acquisition may automatically do this). Designate this spectrum $I_{\text{HKE}}(E)$ (where HKE is 'high kinetic energy').
13. Acquire a survey spectrum between 180 eV and 1500 eV kinetic energy. Designate this spectrum $I_{\text{HKE}}(E)$.
14. Export the $I_{\text{initial}}(E)$, $I_{\text{HKE}}(E)$, and $I_{\text{final}}(E)$ spectra in counts per second and plot them. All three spectra should overlap and have no significant differences in intensity. Also ensure that these spectra have not been intensity calibrated automatically by the acquisition software. Divide

$I_{\text{final}}(E)$ by $I_{\text{initial}}(E)$ and ensure that the ratio is unity across the entire kinetic energy range – if the ratio does not equal unity, then this could indicate an issue with the instrument (sloped ratio), or simply that the X-ray source has not sufficiently equilibrated (offset).

15. Repeat steps 10 to 14 for each of the usual operating modes of the instrument (i.e. pass energies, lens modes, slit / iris settings, etc.).
16. For each operating mode, combine the corresponding $I_{\text{initial}}(E)$, $I_{\text{HKE}}(E)$, and $I_{\text{final}}(E)$ spectra into a single spectrum in counts per second and subtract the average intensity in counts per second above the photon energy D to account for dark noise (or subtract $D(E)$ if variable dark noise was observed in step 7). The combined spectrum $I_{\text{PE}}(E)$ is described mathematically as

$$I_{\text{PE}}(E) = \begin{cases} (I_{\text{initial}}(E) + I_{\text{final}}(E))/2 - D(E), & E < 1195 \text{ eV} \\ I_{\text{HKE}}(E) - D(E), & E \geq 1195 \text{ eV} \end{cases} \quad (4)$$

17. For each operating mode, calculate the LDPE reference spectrum²⁵ $I_{\text{ref}}(E)$ for the relevant instrument geometry $[a, b]$ with the adjustable variable q , which is the kinetic energy offset between $I_{\text{PE}}(E)$ and $I_{\text{ref}}(E)$.
18. For each operating mode divide $I_{\text{PE}}(E)$ by $I_{\text{ref}}(E)$ and adjust q to eliminate sharp features at the C *KLL* edge ($E \sim 280$ eV). The resulting ratio gives $Q_{\text{PE}}(E)$ from which datapoints between 1150 eV and 1220 eV kinetic energy, and values greater than 1440 eV should be removed.
19. Check $Q_{\text{PE}}(E)$ for any issues before continuing, such as intensity steps over the C 1s region, sharp features in the data, uncharacteristic intensity increases near C 1s or C *KLL*, etc. These features could indicate either an issue with the data acquisition (e.g. contaminated LDPE, accidental export of an instrument's old T , etc) or an issue with the spectrometer (e.g. internal scattering, detector saturation, etc). Investigate these issues before continuing, and if necessary, start the calibration procedure again (from step 3).
20. For each operating mode, use power law fits (see S1) in the regions within ~ 100 eV from the ends of $Q_{\text{PE}}(E)$ (i.e. ~ 180 eV to 280 eV and ~ 1340 eV to 1440 eV). Using the power law fits, extrapolate $Q_{\text{PE}}(E)$ below 180 eV (to ~ 100 eV) and above 1440 eV (to ~ 2000 eV). About 10

datapoints above and below the limits of $Q_{PE}(E)$ is enough. This ensures that the subsequent fitting does not drastically diverge beyond the limits of $Q_{PE}(E)$.

21. For each operating mode, fit $Q_{PE}(E)$ with a functional form (equation 2 from Shard & Spencer²²) to obtain the transmission function $T_{PE}(E)$.
22. The transmission function $T_{PE}(E)$ may be used to intensity correct spectra acquired using the same operating parameters: To correct the full spectrum, the intensity at each energy E is divided by $T_{PE}(E)$; this method is necessary if inelastic background analysis is intended. To correct peak areas, the peak area and peak energy, E_P , are evaluated in the raw data and the peak area divided by $T_{PE}(E_P)$ prior to the application of sensitivity factors or other types of analysis.

V. CONCLUSIONS

The results of this inter-laboratory study have shown that the LDPE calibration method can provide an accurate transmission function that agrees with the more established NPL calibration method. Except for two submitted datasets (due to spectrometer issues), every participant's LDPE spectra could be used to determine the function T for their spectrometer; therefore, highlighting the usability and validity of the LDPE intensity calibration method across different laboratories and XPS instruments. We note that the high pass energy datasets result in less systematic deviation compared to the low pass energy datasets, that is $\pm 5.7\%$ for high pass energies and $\pm 8.8\%$ for low pass energies. In short, this is due to a difference in the intensities of the acquired I_{PE} : the high pass energy datasets have a greater signal-to-noise and are less affected by errors in determining the dark noise contribution. In practice, most XPS users will only require an accurate transmission function for their high pass energy acquisition parameters, which will be used to obtain atomic compositions from a survey scan. If the transmission function does not greatly vary over the energy range of the calibration, then the analyser response over a smaller range, such as over a single core level, may be ignored for the purposes of a standalone peak fit. The larger systematic deviation at low pass energies may be reduced by increasing the acquisition time to improve signal-to-noise, although the Q_{PE} will still be highly sensitive to any dark noise

contribution. The inter-laboratory study has also highlighted some advantages and disadvantages to the method in addition to those discussed at the beginning of this article, i.e. *ex-situ* sample preparation without an Ar-ion sputter source, insensitive to carbon contamination and simple mathematically generated reference spectra. A continuous noise-free reference spectrum for LDPE can now be generated for any quartz-monochromated Al K α instrument in the $[a, b]$ geometry. The photoelectron yield from LDPE is low which means that that long acquisition times are required to obtain a suitable signal-to-noise, especially in the HKE region above the C 1s peak. However, this means that LDPE is very unlikely to saturate the detector and cause non-linearity effects in the spectrum even at high X-ray powers, which is an issue that can arise when calibrating with metals like gold, silver, and copper. This means that the operator can calibrate their instrument using their normal acquisition parameters and X-ray source power. The low intensity of the LDPE spectrum means that it is more sensitive to underlying spectrometer issues such as low-level internal scattering or variable dark noise, as these can prevent the fitting of an appropriate T . However, this also means that the LDPE intensity calibration method can potentially be used to diagnose issues with the spectrometer which may be missed when using gold, silver, or copper as the reference material. With these attributes in mind, LDPE appears to be a suitable candidate as a secondary reference material for XPS intensity calibration, and could contribute to the formulation of an international standard as an alternative calibration method.

SUPPLEMENTARY MATERIAL

See supplementary material at [URL will be inserted by AIP Publishing] for: (S1) Protocol for VAMAS inter-laboratory study circulated to participants; (S2) Table showing important experimental information, percentage values of Δ (%) and Σ (%), and comments which describe issues or observations with the participants' datasets or direct feedback from the participants; (S2 continued) Figures showing participants' transmission functions calculated from LDPE and gold; (S3) Supporting figures including: a graphical representation of R_i , Δ and Σ from equations (2) and (3); transmission functions calculated at NPL using LDPE prepared with differing surface topographies.

ACKNOWLEDGEMENTS

We acknowledge funding from the project NMS/ST20 of the National Measurement System by the UK Department for Business, Energy and Industrial Strategy. Thank you to the chair of VAMAS TWA 2, Prof. Ian Gilmore and Dr Charles Clifford for their help initiating the inter-laboratory study. Our sincerest thanks to all the participants of the VAMAS A27 inter-laboratory study; A. P. B. would like to acknowledge that a portion of this research was conducted at the Center for Nanophase Materials Sciences, which is a DOE Office of Science User Facility; D.G.C. gratefully acknowledges support from the US National Institutes of Health grant EB-002027; E. F. S would like to acknowledge the EPSRC grant for the Kratos LiPPS XPS instrument EP/K005138/1.

DATA AVAILABILITY STATEMENT

The data that support the findings of this study are available from the corresponding author upon reasonable request.

REFERENCES

- ¹C. J. Powell, *Micros. Today* **24**, 16-23 (2016).
- ²A. Jablonski and C. J. Powell, *Surf. Sci. Rep.* **47**, 33-91 (2002).
- ³V. I. Nefedov and I. S. Nefedova, *J. Electron Spectrosc.* **107**, 131-137 (2000).
- ⁴J. J. Yeh and I. Lindau, *Atom. Data Nucl. Data* **32**, 1-155 (1985).
- ⁵J. H. Scofield, *J. Electron Spectrosc.* **8**, 129-137 (1976).
- ⁶M. Baker, *Nature* **533**, 353-366 (2016).
- ⁷D. R. Baer and M. H. Engelhard, *J. Surf. Anal.* **26**, 94-95 (2019).
- ⁸M. R. Linford, V. S. Smentkowski, J. T. Grant, C. R. Brundle, P. M. A. Sherwood, M. C. Biesinger, J. Terry, K. Artyushkova, A. Herrera-Gómez, S. Tougaard, W. Skinner, J.-J. Pireaux, C. F. McConville, C. D. Easton, T. R. Gengenbach, G. H. Major, P. Dietrich, A. Thissen, M. Engelhard, C. J. Powell, K. J. Gaskell and D. R. Baer, *Microsc. Microanal.* **26**, 1-2 (2020).
- ⁹International Organization for Standardization, "ISO 15472:2010 — Surface chemical analysis — X-ray photoelectron spectrometers — Calibration of energy scales," Geneva, Switzerland (2010).
- ¹⁰International Organization for Standardization, "ISO 24237:2005 — Surface chemical analysis — X-ray photoelectron spectroscopy — Repeatability and constancy of intensity scale," Geneva, Switzerland (2005).

- ¹¹L. T. Weng, G. Vereecke, M. J. Genet, P. Bertrand and W. E. E. Stone, *Surf. Interface Anal.* **20**, 179-192 (1993).
- ¹²R. Hesse, P. Streubel and R. Szargan, *Surf. Interface Anal.* **37**, 589-607 (2005).
- ¹³R. C. Wicks and N. J. C. Ingle, *Rev. Sci. Instrum.* **80**, 053108 (2009).
- ¹⁴J. Trigueiro, W. Lima, N. Bundaleski and O. M. N. D. Teodoro, *J. Electron Spectrosc.* **222**, 122-132 (2018).
- ¹⁵B. Gruzza, P. Bondot, A. Porte, C. Jardin and G. Gergely, *Acta Phys. Pol. A* **81**, 159-164 (1992).
- ¹⁶C. S. Hemminger, T. A. Land, A. Christie and J. C. Hemminger, *Surf. Interface Anal.* **15**, 323-327 (1990).
- ¹⁷J. Osterwalder, M. Sagurton, P. J. Orders, C. S. Fadley, B. D. Hermsmeier and D. J. Friedman, *J. Electron Spectrosc.* **48**, 55-99 (1989).
- ¹⁸M. Holzweber, A. Lippitz, R. Hesse, R. Denecke, W. S. M. Werner and W. E. S. Unger, *J. Electron Spectrosc.* **233**, 51-56 (2019).
- ¹⁹M. P. Seah, *J. Electron Spectrosc.* **71**, 191-204 (1995).
- ²⁰M. P. Seah, *Surf. Interface Anal.* **20**, 243-266 (1993).
- ²¹M. P. Seah and G. C. Smith, *Surf. Interface Anal.* **15**, 751-766 (1990).
- ²²A. G. Shard and S. J. Spencer, *Surf. Interface Anal.* **51**, 618-626 (2019).
- ²³A. G. Shard, *Surf. Interface Anal.* **46**, 175-185 (2014).
- ²⁴B. P. Reed, S. J. Spencer and A. G. Shard, "VAMAS TWA 2, 2019, sub-project A27 – "intensity calibration for XPS instruments using low-density poly(ethylene)" – protocol for analysis," NPL Report AS 100, National Physical Laboratory, Teddington, UK (2019).
- ²⁵A. G. Shard and B. P. Reed, *unpublished*
- ²⁶A. Herrera-Gomez, *J. Electron Spectrosc.* **182**, 81-83 (2010).
- ²⁷M. P. Seah, I. S. Gilmore and S. J. Spencer, *J. Electron Spectrosc.* **104**, 73-89 (1999).
- ²⁸A. G. Shard, *J. Vac. Sci. Technol. A* **38**, 041201 (2020).

Tables

TABLE I. A list of participants' instruments with corresponding geometry. The terms a , b , and G are detailed in the results section and reference 25. Where angle b is not reported (—), an angle of 180° has been used to calculate G .

Code	Spectrometer Manufacturer / Model	Angle a	Angle b	G (2 s.f.)
AA	Thermo Scientific / Alpha 110	41	180	-0.52
AB	Kratos Analytical / Nova	54.7	180	-0.22
AC	PHI / VersaProbe II	45	180	-0.43
AD	Thermo Scientific / Theta 300	56	45	-0.11
AE	PHI / VersaProbe II	45	180	-0.43
AF	Kratos Analytical / Axis Ultra DLD	60	180	-0.10
AG	Thermo Scientific / Nexsa	60	180	-0.10
AH	PHI / VersaProbe II	45	180	-0.43
AI	Scienta Omicron / UHV System	54.7	—	—
AJ	ULVAC-PHI / Quantum 2000	45	180	-0.43
AK	SPECS / PHOIBOS 150	54	180	-0.23
AL	PHI / VersaProbe II	45	180	-0.43
AM	Kratos Analytical / Axis 165	54	180	-0.23
AN	SPECS / PHOIBOS 150	55	180	-0.21
AO	Kratos Analytical / Axis Ultra DLD	60	180	-0.10
AP	Kratos Analytical / Axis Ultra DLD	60	180	-0.10
AQ	Thermo Scientific / VG ESCALAB 250 Xi	58	180	-0.14
AR	PHI / VersaProbe III	45	180	-0.43
AS	Kratos Analytical / Axis Ultra DLD	60	180	-0.10
AT	SPECS / PHOIBOS 150	55	180	-0.21
AU	PHI / Quantera II Hybrid	45	180	-0.43
AV	Surf Sci / SSX-100	75.5	120	0.37
AW	Kratos Analytical / Axis Supra DLD	60	180	-0.10
AX	SPECS / PHOIBOS 150	55	180	-0.21
AY	ULVAC-PHI / Quantera SXM	45	180	-0.43
AZ	Thermo Scientific / Nexsa	60	180	-0.10
AZ	Thermo Scientific / VG ESCALab 250Xi	58	180	-0.14
BA	Kratos Analytical / Axis Supra	60	180	-0.10
BB	PHI / Quantera Hybrid	45	—	—
BC	ULVAC-PHI / Quantera SXM	45	180	-0.43
BD	SSI / S-Probe	55	—	—
BE	Kratos Analytical / Axis Ultra DLD	60	180	-0.10
BF	Thermo Scientific / ESCALAB Xi+	58	180	-0.14
BG	Thermo Scientific / VG ESCALAB 220i-XL	54.7	180	-0.22
BG	ULVAC-PHI / ESCA 5800	60	54.7	0.01
BH	Kratos Analytical / Axis Ultra DLD	60	180	-0.10
BI	Kratos Analytical / Axis Ultra DLD	60	180	-0.10

This is the author's peer reviewed, accepted manuscript. However, the online version of record will be different from this version once it has been copyedited and typeset.
PLEASE CITE THIS ARTICLE AS DOI: 10.1116/6.0000577

TABLE II. Average values of Δ (%) and Σ (%) for all participants' datasets, low pass energy datasets, and high pass energy datasets.

Dataset Grouping	Δ (%)	Σ (%)
All participants	3.0 %	\pm 6.5 %
Low pass energy (< 90 eV)	$-$ 4.9 %	\pm 8.8 %
High pass energy (\geq 90 eV)	5.8 %	\pm 5.7 %

Figure Captions

FIG. 1. Geometry correction to the LDPE I_{ref} . In panels (A) and (B), the black dashed line denotes an axis break in kinetic energy where the low KE region scale is on the left, and the high KE region scale is on the right. (A) The two calibrated LDPE spectra from the VAMAS participants data. The $[45^\circ, 180^\circ]$ spectrum is shown in black, and the $[60^\circ, 180^\circ]$ spectrum is shown in red. (B) Geometry-corrected LDPE I_{ref} calculated using the method described by Shard & Reed.²⁵ A reference spectrum for any combination of a and b can be calculated. (C) The effect of using the new geometry correction on Q_{PE} for participant AA. The red dots show the values of Q_{PE} when calculated using a LDPE I_{ref} for the $[60^\circ, 180^\circ]$ geometry, and the black dots show the values of Q_{PE} when the correct $[41^\circ, 180^\circ]$ geometry LDPE I_{ref} is used.

FIG. 2. Transmission function correction using LDPE and Au spectra from participants AX (red), AJ (blue), AI (green), and AT (orange). (A, B) LDPE and Au spectra normalised to the C KLL and Au NNV Auger peaks respectively. (C, D) Q (dots) and T (lines) calculated using LDPE and Au intensity calibration methods respectively. (E, F) Intensity-corrected LDPE and Au spectra using T_{PE} and T_{Au} respectively.

FIG. 3. The values of Δ (%) and $\pm \Sigma$ (%) for all the participants' datasets, where the red diamonds (◆) refer to high pass energy datasets (≥ 90 eV), and the blue squares (■) refer to low pass energy data sets (< 90 eV). The dashed lines show the average values of Δ (%) and $\pm \Sigma$ (%) from all datasets (black), high pass energies (red), and low pass energies (blue). The pink dots (●) refer to datasets which have been replaced with newer data from the same participant, and do not contribute to the average values of Δ (%) and $\pm \Sigma$ (%).

FIG. 4. The Q (dots) and T (lines) calculated from LDPE (green) and Au (orange) respectively for participants AN-0, BB, and BF-7, which have been highlighted in Figure 3.

FIG. 5. (A) The variable dark noise $D(E)$ recorded by participant AX (black dots) fitted with an arbitrary function $f(E)$ (red line). (B) The effect of correcting $D(E)$ shown by participant AX. The difference in Q_{PE} when


corrected with a constant dark noise value (black +) or $f(E)$ (red ×). In the high kinetic energy, low count rate region above the C 1s, the difference in Q_{PE} is ~4.9 %.

FIG. 6. Q_{PE} and Q_{Au} from participant BD showing significant internal scattering and intensity loss in the low kinetic energy region.


FIG. 7. Q_{PE} and Q_{Au} from participant AD-2 demonstrating the significant effect that contamination on the LDPE produces. Other elements (i.e. O, Si, and S) on the LDPE cause steps, discontinuities, and deviations in Q_{PE} .

FIG. 8. **(A)** Q_{Au} (red dots) and T_{Au} (black line) from participant AF. The marked positions show an uncharacteristic increase in the spectrum intensity near the gold peaks which may affect the fitting of T_{Au} . **(B)** Plotting the actual I_{Au} from participant AF (red dots) against the expected I_{Au-ref} (intensity corrected using T_{PE} , and corrected for average offset) (black line) reveals detector non-linearity above 300 kcps.


This is the author's peer reviewed, accepted manuscript. However, the online version of record will be different from this version once it has been copyedited and typeset.
PLEASE CITE THIS ARTICLE AS DOI: 10.1116/6.0000577

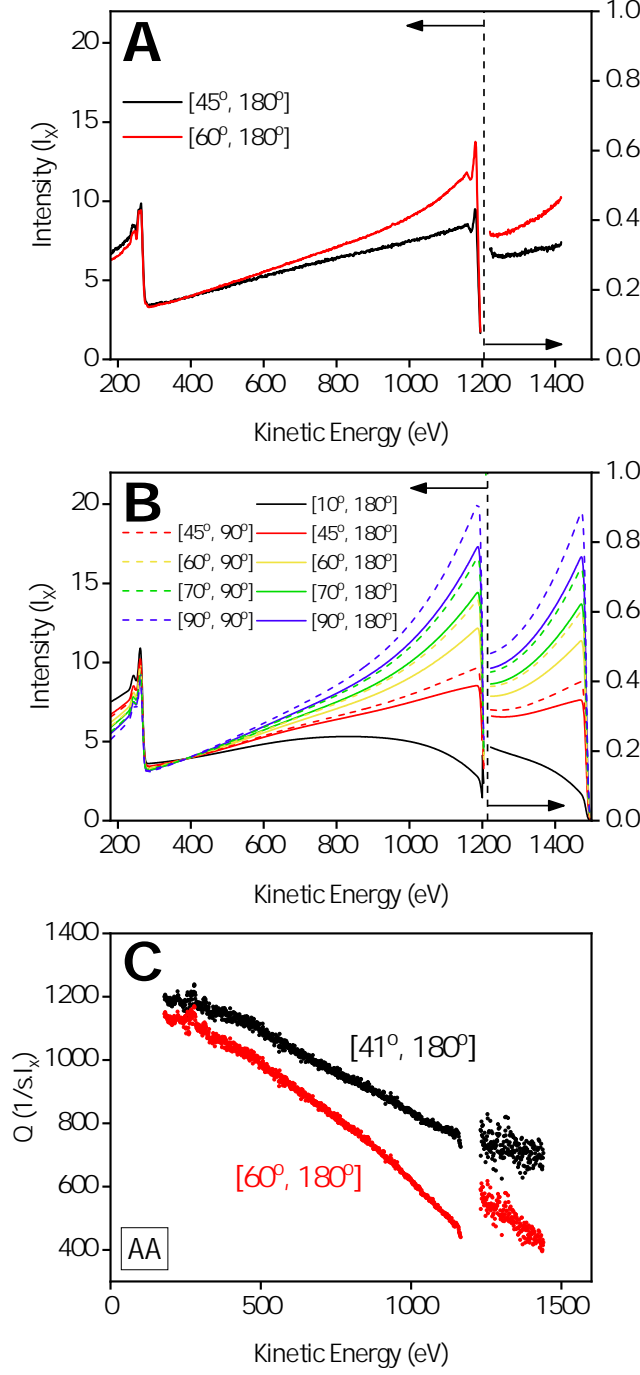


This is the author's peer reviewed, accepted manuscript. However, the online version of record will be different from this version once it has been copyedited and typeset.
PLEASE CITE THIS ARTICLE AS DOI: 10.1116/1.5000057

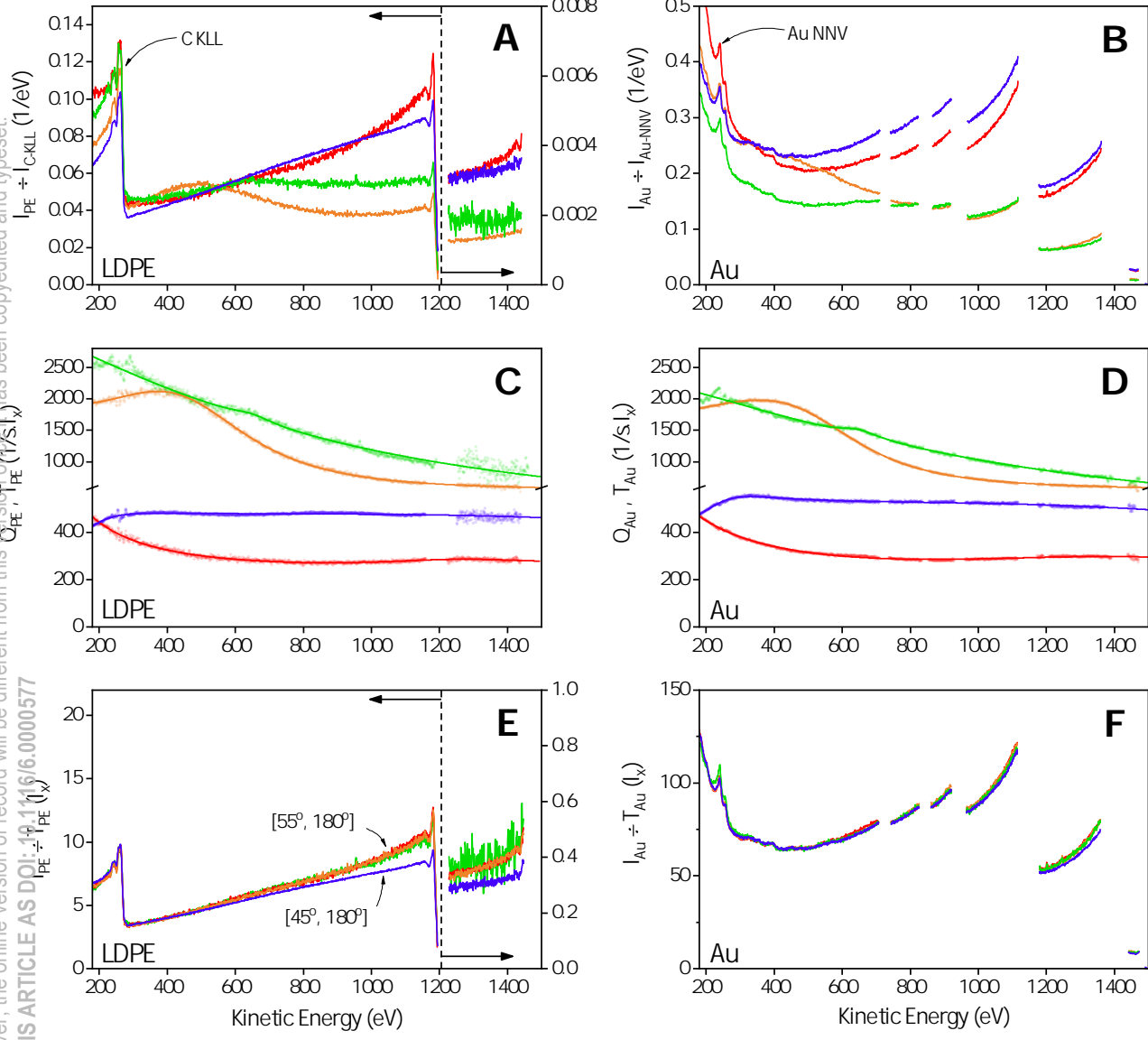


This is the author's peer reviewed, accepted manuscript. However, the online version of record will be different from this version once it has been copyedited and typeset.
PLEASE CITE THIS ARTICLE AS DOI: 10.1116/6.0000577

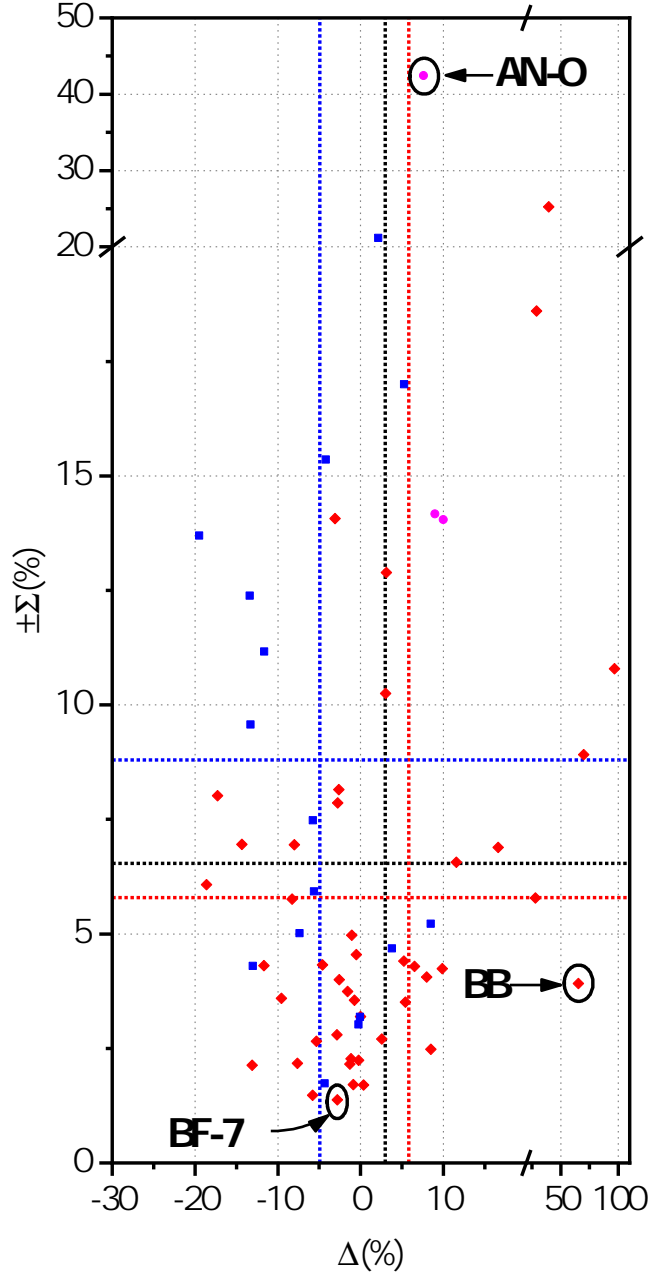


This is the author's peer reviewed, accepted manuscript. However, the online version of record will be different from this version once it has been copyedited and typeset.
PLEASE CITE THIS ARTICLE AS DOI: 10.1116/6.0000577

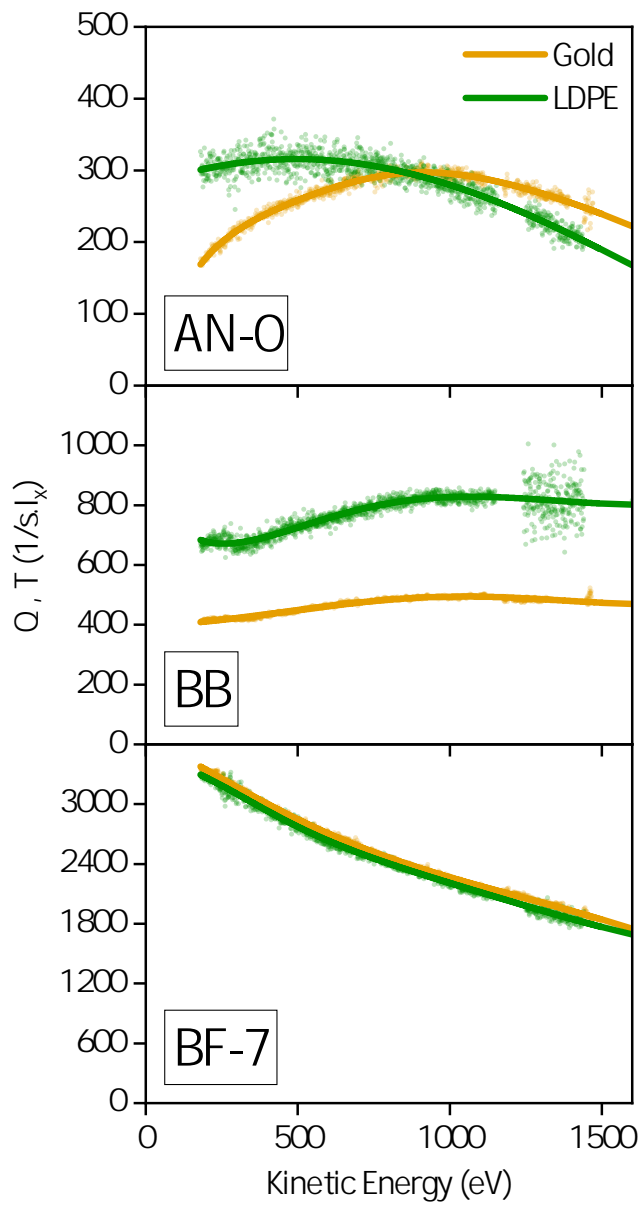


This is the author's peer reviewed, accepted manuscript. However, the online version of record will be different from this version once it has been copyedited and typeset.
PLEASE CITE THIS ARTICLE AS DOI: 10.1116/6.0000577

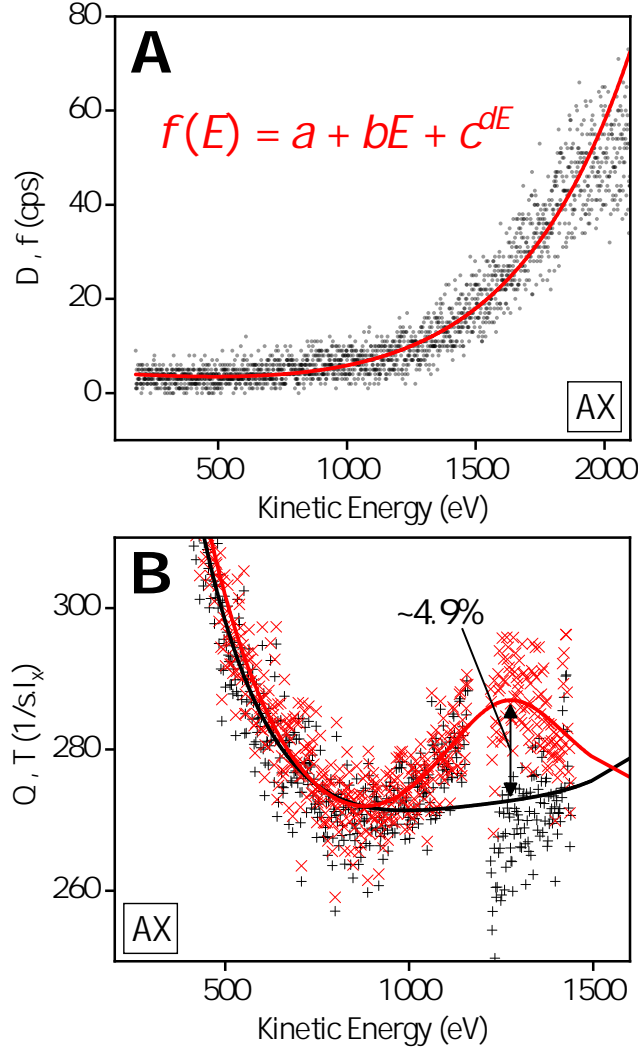

This is the author's peer reviewed, accepted manuscript. However, the online version of record will be different from this version once it has been copyedited and typeset.
PLEASE CITE THIS ARTICLE AS DOI: 10.1116/6.0000577

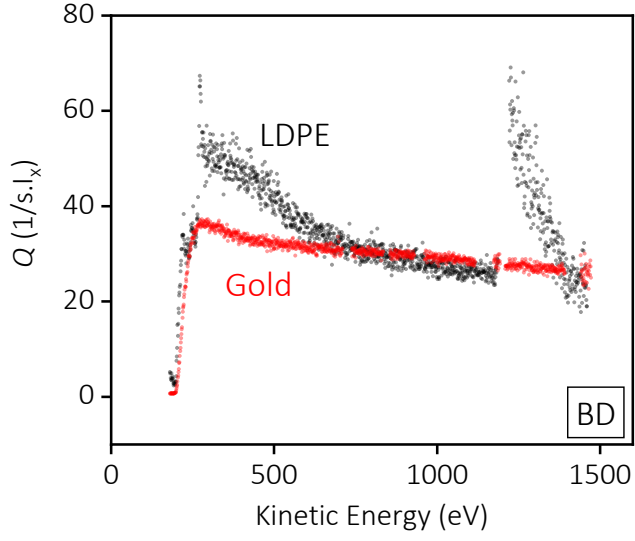

This is the author's peer reviewed, accepted manuscript. However, the online version of record will be different from this version once it has been copyedited and typeset.
PLEASE CITE THIS ARTICLE AS DOI: 10.1116/6.0000577

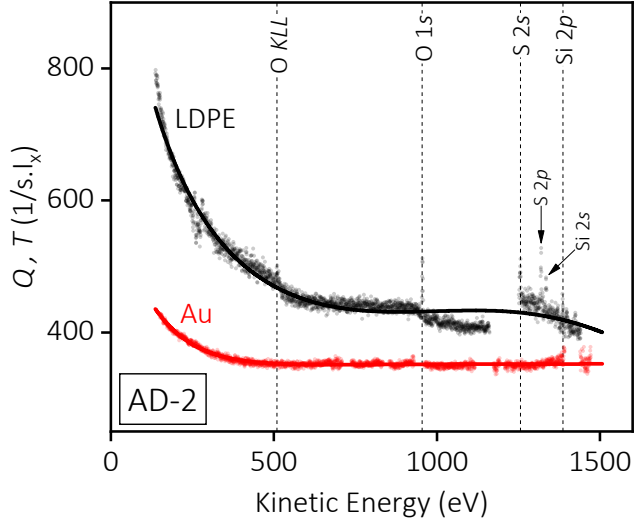

This is the author's peer reviewed, accepted manuscript. However, the online version of record will be different from this version once it has been copyedited and typeset.
PLEASE CITE THIS ARTICLE AS DOI: 10.1116/6.0000577

